

Datasheet for 200-301-NF4**ModDetect® 2'-O-Methoxyethyl (2'MOE) Clone MOEC****Overview**

Description:	ModDetect® 2'-O-Methoxyethyl (2'MOE) Clone MOEC - 200-301-NF4
Item No.:	200-301-NF4
Size:	100 µg
Applications:	ELISA, IF, Multiplex
Host Species:	Mouse

Product Details

Background:	Current studies focusing on RNA therapeutics use as potential drugs antisense oligonucleotides (ASO), short interfering RNA (siRNA), and micro RNA (miRNA). The choice of oligonucleotide chosen for a specific treatment is dependent on many factors. However, all oligo therapeutics must be modified in some manner, as unmodified nucleotides are difficult to identify and are highly susceptible to degradation by endogenous nucleases. Modifications can be created within the sugar-phosphate backbone (first generation), the ribose sugar (second and third generation), or at multiple locations (third generation). Antibodies to modified oligonucleotides are useful for basic immunoassays including, ELISA, IHC, and IF. In addition, several types of toxicology assays (immunogenicity assays), such as antibody-drug assays (ADA) or other pharmacokinetic (PK) and pharmacodynamic (PD) studies require oligonucleotide antibodies as analytical tools for general research and preclinical trial experiments. Antibodies to an oligonucleotide therapeutic allow researchers to assess cellular uptake, tissue distribution, and can serve as a positive control for immunogenicity studies.
Synonyms:	2'-O-Methoxyethyl, 2'MOE, ASO, anti-sense oligonucleotide, siRNA, RNA therapeutic, gapmer, aptamer
Host Species:	Mouse
Clonality:	Monoclonal
Clone ID:	MOEC
Format:	IgG1

Target Details

Immunogen Type:	Other
------------------------	-------

Immunogen:	ModDetect® 2'-O-Methoxyethyl (2'MOE) Clone MOEC was prepared from cell culture supernatant produced by repeated immunizations with a proprietary oligo sequence.
Purity/Specificity:	This Protein A purified ModDetect® 2'-O-Methoxyethyl (2'MOE) Clone MOEC is directed against 2'MOE modified ribose sugar. This product was Protein A purified from monospecific supernatant.

Application Details

Tested Applications:	ELISA, IF, Multiplex
Application Note:	ModDetect® 2'-O-Methoxyethyl (2'MOE) Clone MOEC has been validated for use in ELISA. In antigen-down or immunometric ELISAs, the ModDetect® 2'-O-Methoxyethyl (2'MOE) Clone MOEC preferentially detects 2'MOE modified ribose sugar groups over native deoxyribose or ribose sugars. It is also likely to be functional in IF and IHC assays, however, its sensitivity, specificity, and conditions for optimal use in these assays has not yet been determined. All assays should be optimized by the user.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
ELISA:	1:1000

Formulation

Physical State:	Liquid (sterile filtered)
Concentration:	3.001 by UV absorbance at 280 nm
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2
Preservative:	0.01% (w/v) Sodium Azide
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free

Shipping & Handling


Shipping Condition:	Dry Ice
Storage Condition:	Store vial at -20° C or below prior to opening. Store the vial at -20°C or below after dilution. Avoid cycles of freezing and thawing.
Expiration:	Expiration date is one (1) year from date of receipt.

Images



Vial


ModDetect® 2'-O-Methoxyethyl (2'MOE) Clone MOEC

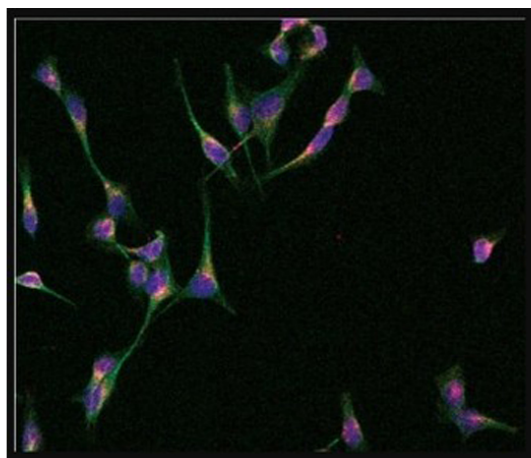

Immunofluorescence Microscopy

Colocalization of Cy5-conjugated PS-MOE ASO with ModDetect® Anti-MOE antibody.

HeLa cells were cultured and treated with 1 μ M PS-MOE-modified ASO conjugated to the fluorochrome Cy5 by gymnotic uptake. After fixation with paraformaldehyde, cells were stained with DAPI (blue), anti-alpha tubulin (green), and ModDetect® anti-MOE antibody clone MOE-C (red) [p/n 200-301-NF4]. The fluorescence from Cy5-conjugated ASO (magenta) is shown to colocalize with anti-MOE antibody staining in the merged image.

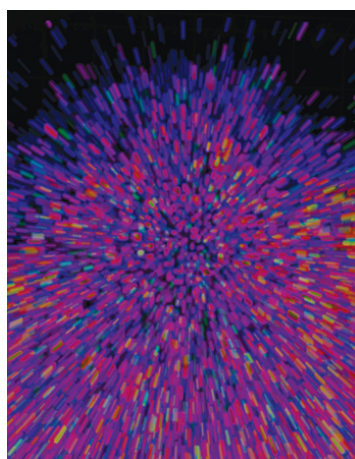

Immunocytochemistry

Immunocytochemistry using ModDetect®. This image, featured on the cover of OTS Nucleic Acid Therapeutics, Vol. 35, No. 4, August 2025, features a three-dimensional reconstruction with pseudo-color rendering of a spheroid culture with staining for cell nuclei and immunocytochemistry using antibodies to detect endosomes and the phosphorothioate (PS) modification of an antisense oligonucleotide. These cultures were used as part of a benchmarking study in cells and tissue to assess the ModDetect® series of antibodies raised against PS and 2'-MOE.

Image provided by Peter Oliver of the Medical Research Council Nucleic Acid Therapy Accelerator, Research Complex at Harwell, Didcot, UK.

References

- Chimento DP et al. Bioanalytical Assays for Oligonucleotide Therapeutics: Adding Antibody-Based Immunoassays to the Toolbox as an Orthogonal Approach to LC-MS/MS and Ligand Binding Assays. *Nucleic Acid Ther.* (2025)
- Fial I et al. Characterizing Antibodies Targeting Antisense Oligonucleotide Phosphorothioate and 2'-O-Methoxyethyl Modifications for Intracellular Trafficking and Biodistribution Studies. *Nucleic Acid Ther.* (2025)

Disclaimer

THE MATERIALS PROVIDED WITH THIS CERTIFICATE OF ANALYSIS (MATERIALS) ARE FOR RESEARCH USE ONLY AND SUBJECT TO A MATERIAL EVALUATION AGREEMENT.

BY USING THIS PRODUCT, THE USER AGREES TO BE BOUND BY THE TERMS OF THE MATERIAL EVALUATION AGREEMENT. If the user is not willing to accept the terms of the agreement and the product is unused, Rockland will accept return of the unused product and provide the user with a full refund. Pursuant to the agreement, the user may use the product for research use only.

NO TRANSFER OR COMMERCIAL USE OF THIS PRODUCT IS ALLOWED. Commercial use means a use of the product in exchange for consideration, including use in the manufacture of a product, in a clinical trial or other testing regulated by a government agency (e.g. FDA, EMEA, EPA, etc.) including ADA and clearance assays, use of the product in the provision of services, and the resale of the product. Rockland reserves all other rights not permitted under the license granted herein, including the rights which are embodied in other products or materials derived from the use of the Materials or know-how embodied in the Materials.

NO RIGHT IS GRANTED UNDER THE LICENSE TO MODIFY OR OTHERWISE CREATE DERIVATIVES OR VARIATIONS OF THE PRODUCT OR TO USE ROCKLAND'S KNOW-HOW OR OTHER INTELLECTUAL PROPERTY FOR OTHER USES. No other use of this product is authorized without the express written prior consent of Rockland, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, therapeutic uses, or any type of consumption by or application to humans or animals.

For uses outside the agreement including any diagnostic, therapeutic or commercial use, please contact Rockland for supply and licensing information at licensing@rockland.com.