

Datasheet for 600-401-105S**Collagen Type III Antibody****Overview**

Description:	Anti-Collagen Type III (RABBIT) Antibody - 600-401-105S
Item No.:	600-401-105S
Size:	25 µL
Applications:	Dot Blot, IHC, WB, ELISA, FC, IF, Multiplex
Reactivity:	Human, Bovine, Pig
Host Species:	Rabbit

Product Details

Background:	Rockland produces highly active antibodies and conjugates to collagens. Collagens are highly conserved throughout evolution and are characterized by an uninterrupted "Glycine-X-Y" triplet repeat that is a necessary part of the triple helical structure. For these reasons, it is often extremely difficult to generate antibodies with specificities to collagens. The development of 'type' specific antibodies is dependent on NON-DENATURED three-dimensional epitopes. Rockland extensively purifies collagens for immunization from human and bovine placenta and cartilage by limited pepsin digestion and selective salt precipitation. This preparation results in a native conformation of the protein. Antibodies are isolated from rabbit antiserum and are extensively cross-adsorbed by immunoaffinity purification to produce 'type' specific antibodies. Greatly diminished reactivity and selectivity of these antibodies will result if denaturing and reducing conditions are used for SDS-PAGE and immunoblotting. Ideal for investigators involved in Cell Biology, Signal Transduction and Stem Cell research.
Synonyms:	rabbit anti-Collagen Type III antibody, Collagen type III alpha 1 antibody, Collagen type III alpha antibody, EDS4A antibody, Ehlers Danlos syndrome type IV, autosomal dominant antibody, Fetal collagen antibody, COL3A1, Collagen alpha-1 (III) chain
Host Species:	Rabbit
Clonality:	Polyclonal
Format:	IgG

Target Details

Gene Name:	COL3A1
-------------------	--------

Reactivity:	Human, Bovine, Pig
Immunogen Type:	Native Protein
Immunogen:	Collagen Type III from human and bovine placenta
Purity/Specificity:	Collagen III Antibody has been prepared by immunoaffinity chromatography using immobilized antigens. Some class-specific anti-collagens may be specific for three-dimensional epitopes which may result in diminished reactivity with denatured collagen or formalin-fixed, paraffin embedded tissues. This antibody reacts with most mammalian Type III collagens and has expected cross-reactivity with Type I and negligible cross reactivity with Type II, IV, V or VI collagens. Non-specific cross-reaction of anti-collagen antibodies with other human serum proteins or non-collagen extracellular matrix proteins has not been tested.
Relevant Links:	Anti-Collagen IHC Protocol UniProtKB - P02461 NCBI - NP_000081.1 GeneID - 1281

Application Details

Tested Applications:	Dot Blot, IHC, WB
Suggested Applications:	ELISA, FC, IF, Multiplex (Based on references)
Application Note:	Anti-Collagen Type III has been tested by dot Blot, western blot, and IHC and is useful for indirect trapping ELISA for quantitation of antigen in serum using a standard curve, immunoprecipitation, native (non-denaturing, non-dissociating) PAGE, immunohistochemistry, and western blotting for highly sensitive qualitative analysis.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
ELISA:	1:5,000 - 1:50,000
FC:	User Optimized
IF:	User Optimized
IHC:	1:50 - 1:200
IP:	1:100
WB:	1:1,000 - 1:10,000

Formulation


Physical State:	Liquid (sterile filtered)
------------------------	---------------------------

Concentration:	1.16 mg/ml by UV absorbance at 280 nm
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2
Preservative:	0.01% (w/v) Sodium Azide
Stabilizer:	None

Shipping & Handling


Shipping Condition:	Dry Ice
Storage Condition:	Store vial at -20° C or below prior to opening. This vial contains a relatively low volume of reagent (25 µL). To minimize loss of volume dilute 1:10 by adding 225 µL of the buffer stated above directly to the vial. Recap, mix thoroughly and briefly centrifuge to collect the volume at the bottom of the vial. Use this intermediate dilution when calculating final dilutions as recommended below. Store the vial at -20°C or below after dilution. Avoid cycles of freezing and thawing.
Expiration:	Expiration date is one (1) year from date of receipt.

Images

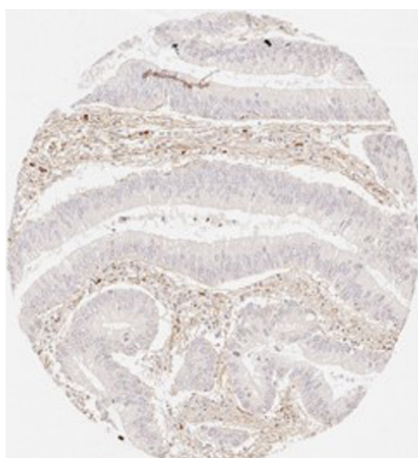


Immunohistochemistry

Immunohistochemistry results of Rabbit Anti-Collagen Type I Antibody. Tissue: human stroma in colorectal cancer. Fixation: FFPE. Antigen Retrieval: HIER using Tris-EDTA-citrate buffer pH 7.8 for 5 min. Blocking: Peroxidase-Blocking Solution for 10 min. Primary Antibody: Anti-Collagen Type I (p/n 600-401-103-0.1) at 1:15 for 1 hr at 37 °C. Secondary Antibody: Dako REAL EnVision Detection Kit, Polymer-HRP, Rabbit/Mouse. Counterstain: Hematoxylin for 15 sec. Substrate: DAB-Chromogen, Rabbit/Mouse. Staining/Results: Fibrillar collagen III staining of the stroma in a colorectal cancer. Independently Validated by antibodies-online GmbH (p/n ABIN7565873/ ABIN5596830/ ABIN5596829) courtesy of MS Validated Antibodies.

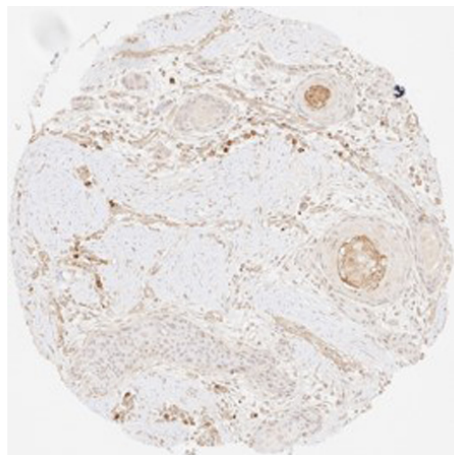


Immunohistochemistry

Immunohistochemistry results of Rabbit Anti-Collagen Type III Antibody. Tissue: human oral cavity. Fixation: FFPE. Antigen Retrieval: HIER using Tris-EDTA-citrate buffer pH 7.8 for 5 min. Blocking: Peroxidase-Blocking Solution for 10 min. Primary Antibody: Anti-Collagen Type III (p/n 600-401-105-0.1) at 1:15 for 1 hr at 37 °C. Secondary Antibody: Dako REAL EnVision Detection Kit, Polymer-HRP, Rabbit/Mouse. Counterstain: Hematoxylin for 15 sec. Substrate: DAB-Chromogen, Rabbit/Mouse. Staining/Results: Distinct fibrillar collagen III staining of the stroma in squamous cell carcinoma of the oral cavity. Independently Validated by antibodies-online GmbH (p/n ABIN7565873/ ABIN5596830/ ABIN5596829) courtesy of MS Validated Antibodies.

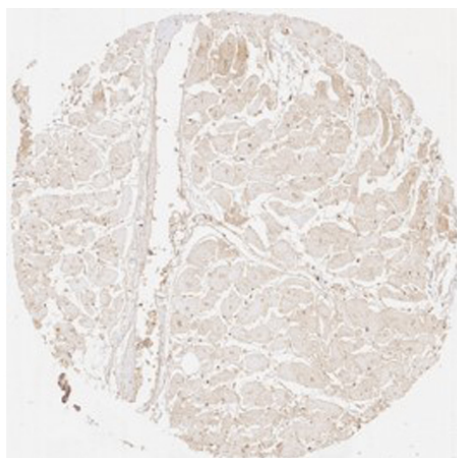


Immunohistochemistry

Immunohistochemistry results of Rabbit Anti-Collagen Type III Antibody. Tissue: human heart muscle. Fixation: FFPE. Antigen Retrieval: HIER using Tris-EDTA-citrate buffer pH 7.8 for 5 min. Blocking: Peroxidase-Blocking Solution for 10 min. Primary Antibody: Anti-Collagen Type III (p/n 600-401-105-0.1) at 1:15 for 1 hr at 37 °C. Secondary Antibody: Dako REAL EnVision Detection Kit, Polymer-HRP, Rabbit/Mouse. Counterstain: Hematoxylin for 15 sec. Substrate: DAB-Chromogen, Rabbit/Mouse. Staining/Results: Distinct fibrillar collagen III staining surrounding each heart muscle cell. Independently Validated by antibodies-online GmbH (p/n ABIN7565873/ABIN5596830/ABIN5596829) courtesy of MS Validated Antibodies.

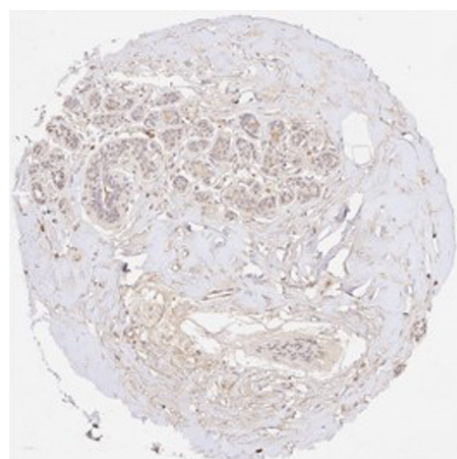


Immunohistochemistry

Immunohistochemistry results of Rabbit Anti-Collagen Type III Antibody. Tissue: human non-cancerous breast tissue. Fixation: FFPE. Antigen Retrieval: HIER using Tris-EDTA-citrate buffer pH 7.8 for 5 min. Blocking: Peroxidase-Blocking Solution for 10 min. Primary Antibody: Anti-Collagen Type III (p/n 600-401-105-0.1) at 1:225 for 1 hr at 37 °C. Secondary Antibody: Dako REAL EnVision Detection Kit, Polymer-HRP, Rabbit/Mouse. Counterstain: Hematoxylin for 15 sec. Substrate: DAB-Chromogen, Rabbit/Mouse. Staining/Results: Fibrillar collagen III staining in non-cancerous breast tissue showing considerable sclerosis. Independently Validated by antibodies-online GmbH (p/n ABIN7565873/ ABIN5596830/ ABIN5596829) courtesy of MS Validated Antibodies.

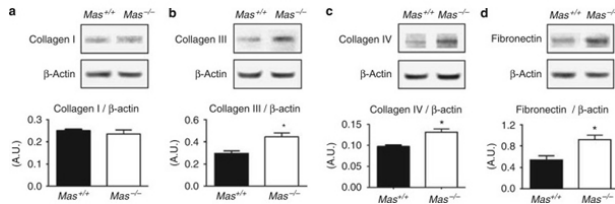


Western Blot

Immunoblotting of extracellular matrix (ECM) proteins in kidneys of Mas^{+/+} and Mas^{-/-} animals. (a) Immunoblotting shows no difference of Collagen I (p/n 600-401-103) expression in Mas^{+/+} and Mas^{-/-} mice kidneys. Significant increases in (b) Collagen III (p/n 600-401-105), (c) Collagen IV (p/n 600-401-106), and (d) fibronectin (p/n 600-401-117) expression were detected by comparing immunoblots of Mas^{-/-} mouse kidneys with those of Mas^{+/+} controls. Each band represents one mouse kidney from either Mas^{+/+} or Mas^{-/-} mice. Data are shown as the mean±s.e.m. *P<0.05. A.U. indicates arbitrary unit. Fig 5. PMID: 19262461

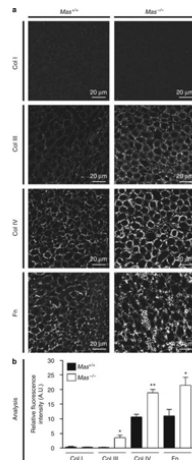


Immunofluorescence Microscopy

Immunofluorescence of extracellular matrix (ECM) proteins in the medulla of kidneys from Mas^{+/+} (left column) and Mas^{-/-} (right column) mice. (a) Fluorescence (Cy3-labeled anti-rabbit IgG) reveals the immunolabeling of ECM proteins. Expression of type III collagen (Col III) (p/n 600-401-105), type IV collagen (Col IV) (p/n 600-401-106), and fibronectin (Fn) (p/n 600-401-117) were increased in the medulla of Mas^{-/-} compared with that of Mas^{+/+} mice, whereas the expression of type I collagen (Col I) (p/n 600-401-103) was unaltered. (b) Quantification of ECM proteins in the medulla of Mas^{+/+} and Mas^{-/-} mice. Data are shown as mean±s.e.m. *P<0.05; **P<0.01. A.U. indicates arbitrary unit. Fig 3. PMID: 19262461

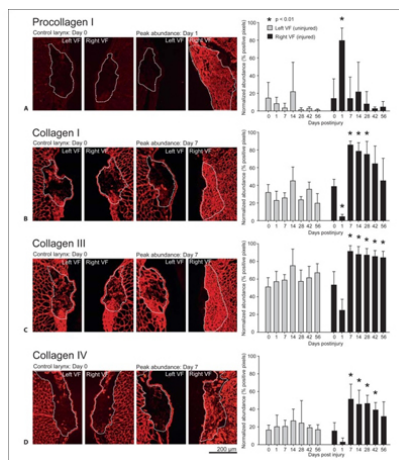

Immunohistochemistry

Representative IHC images and quantitative analysis of procollagen type I (A), collagen type I (B), collagen type III (C), and collagen type IV (D) in uninjured control and injured vocal fold LP over time. Normalized abundance is expressed with respect to total LP area (indicated by a white dashed line). Error bars represent standard deviations and asterisks denote statistically significant differences (p < 0.01) from the control group (day 0). VF = Vocal fold. Scale bar: 200 μm. Fig 2. PMID: 20516667

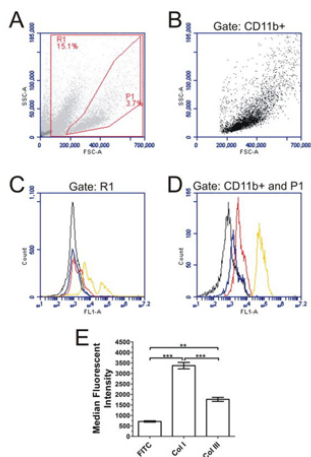



Flow Cytometry

ACK-treated spleen cells were cultured for 5 days at 3.5×10^5 cells per well in the presence of IL-13 and M-CSF. Adherent cells were removed by trypsin-EDTA treatment and stained for the presence of CD11b. The cells were then fixed, permeabilized, stained with rabbit polyclonal antibodies, and analyzed by flow cytometry. A) Forward and side scatter characteristics of 5 day cultured spleen cells. (B) Forward and side scatter analysis of CD11b+ cells. C and D) Histograms show fluorescence intensity of FITC-conjugated goat F(ab')₂ anti-rabbit 2° (black line) compared to collagen I (red line), collagen III (blue line), and as a positive control syk (yellow line). Flow cytometry plots are representative of 3 independent experiments. C) When gating the entire live cell population R1 (from A), collagen was nearly undetectable, but when gating (D) for CD11b+ cells and region P1 (from A and B) there was an increase in the levels of both collagen I and III. E) Compared to FITC control, CD11b+ cells had a significant increase in median fluorescent intensity for collagen I and III staining (1-way ANOVA, Dunnett's test). **, $p < 0.01$; ***, $p < 0.001$. Fig 2. PMID: 20888336

Western Blot

B: representative Western blot of fibronectin (FN) and collagen type III (Col III) levels in C3H 10T1/2 fibroblasts treated with TGF- β 1 5 ng/ml. Tubulin was used as the loading control. Fig 3. PMID: 29351412

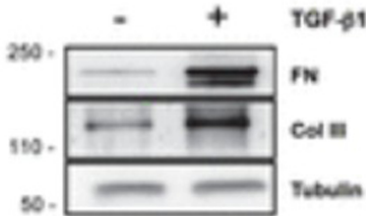B

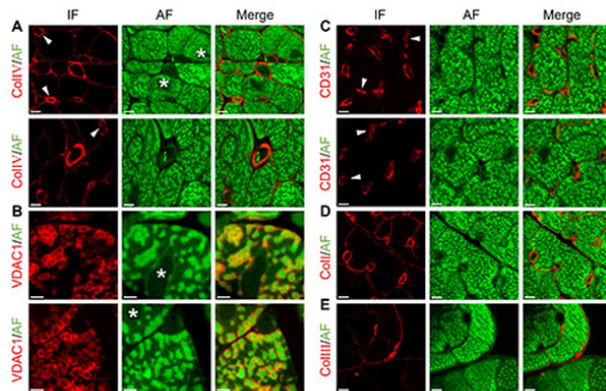


Immunofluorescence Microscopy

Major autofluorescent structures in formaldehyde-fixed mouse myocardial samples identified by IF stainings of (A) Col IV, (B) VDAC1, (C) CD31, (D) Col I (p/n 600-401-103, and (E) Col III (p/n 600-401-105). Images were acquired at cardiomyocyte cross-sectional areas. IF staining (red) allows the localization of basement membrane (Col IV), mitochondria (VDAC1), microvascular endothelial cells (CD31), and interstitial collagen network (Col I and Col III), facilitating the identification of observable structures in AF images (green) acquired at CH1 (525±25 nm). Oval-to-rectangular nuclei (star signs) are observed within cardiomyocytes, and erythrocytes are present within some capillaries (arrowheads). Scale bars: A,C-E) 5 µm; B) 2 µm.

Fig 3.


PMID: 37781779

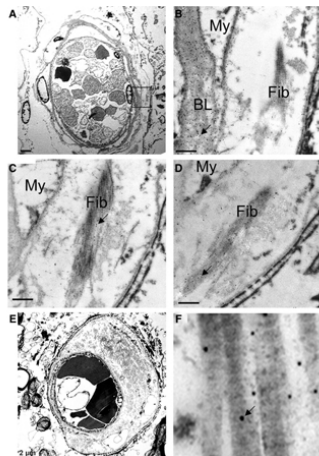

Immunofluorescence Microscopy


Immuno-electron microscopy of collagens I, III, and IV in small vessel disease (SVD). Transmission electron micrographs of small arteries in white matter from older persons with SVD. A–D, Semi-serial sections of a small arteriole, shown at low magnification in A and at higher magnification in B–D with immunogold labeling for collagen-IV (p/n 600-401-106) (B), collagen-I (p/n 600-401-103) (C), or collagen-III (p/n 600-401-105) (D). The region shown at higher magnification in B–D is marked in A (box). Arrows show examples of immunogold particles. In B, collagen-IV labeling shows numerous gold particles over the homogenous basal laminae (BL), behind the endothelium. Mural collagen fibrils (Fib) and myocytes (My) are unlabeled. C, Collagen-I, in contrast to collagen-IV, is mainly localized to the collagen fibrils. D, Collagen-III also labels the banded collagen fibrils. A few gold particles are found over the basal laminae. E and F, Another small artery with severe fibrosis. Low magnification view (E) shows heavy deposits of fibrillar collagens in the vessel wall. In higher magnification, these are heavily labeled for collagen-I (F). Scale bars: 2 µm (A and E), 0.5 µm (B–D), and 0.2 µm (F). Fig 3.

PMID: 36205142

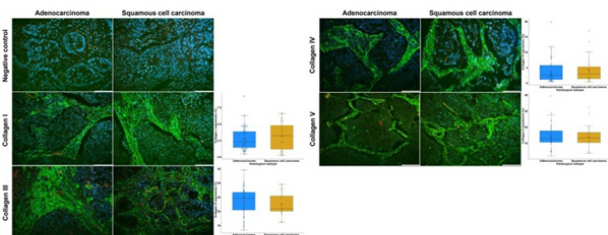


Immunofluorescence Microscopy

Immunofluorescence expression of collagens in the extracellular matrix in lung adenocarcinoma and lung squamous cell carcinoma. The panel shows the negative control, collagen I, collagen III, collagen IV and collagen V. Also as presented the respectively box-plots show the differences of expression between the histological subtypes. The red square in the graphs represents the average expression of the marker for each of the groups. Magnification: 400×. Scale: 100 μm. Fig. 2. PMID: 39938232

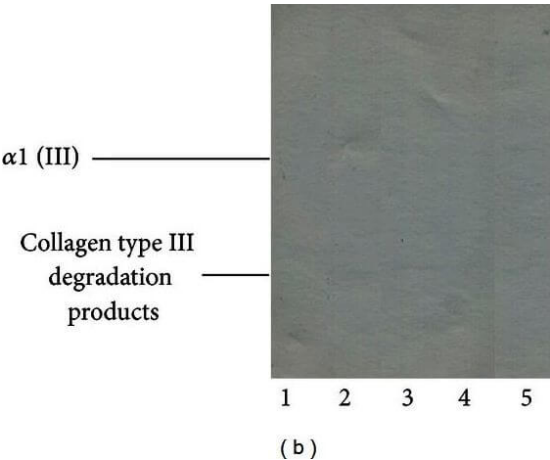


Western Blot

(a) Typical electropherogram of α (I) subunits, oligomers β/γ (I), and collagen type I degradation products, released from healthy and burnt skin. Collagenous components released from tissue samples using pepsin were submitted to 4–15% gradient SDS-PAGE, in nonreducing conditions. (b) Interferences of collagen type III, in electrophoretic profiles of collagen type I components extracted from healthy and burned skin. Collagen components were submitted to electrophoresis in the absence of dithiothreitol (reducing disulfide bonds) and subsequently—after electrotransfer to Immobilon—subjected to reaction with collagen type III antibodies. Lane 1: components of collagen type I isolated from healthy skin. Lane 2: components of collagen type I isolated from burned skin treated with propolis. Lane 3: components of collagen type I isolated from burned skin treated with AgSD. Lane 4: components of collagen type I isolated from burned skin treated with propolis vehicle. Lane 5: components of collagen type I isolated from burned skin treated with NaCl. Figure provided by CiteAb. Source: Evid Based Complement Alternat Med, PMID: 23781260.

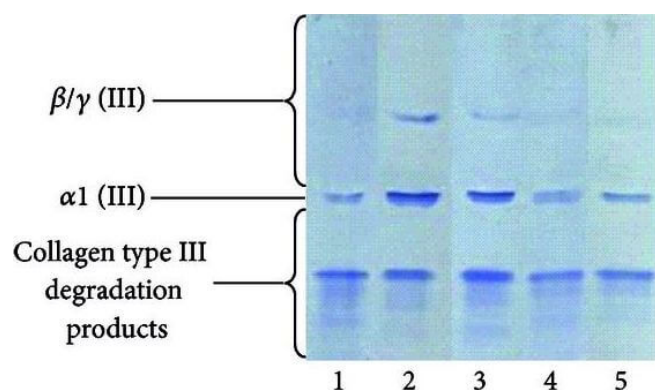


Western Blot

Typical electropherogram of α (III) subunits, oligomers β/γ (III), and collagen type III degradation products, released from healthy and burned skin. Collagenous components released from tissue samples using pepsin were submitted to 4–15% gradient SDS-PAGE, western blotted, and probed with collagen type III antibodies. Lane 1: components of collagen type I isolated from healthy skin. Lane 2: components of collagen type I isolated from burned skin treated with propolis. Lane 3: components of collagen type I isolated from burned skin treated with AgSD. Lane 4: components of collagen type I isolated from burned skin treated with propolis vehicle. Lane 5: components of collagen type I isolated from burned skin treated with NaCl. Figure provided by CiteAb. Source: Evid Based Complement Alternat Med, PMID: 23781260.

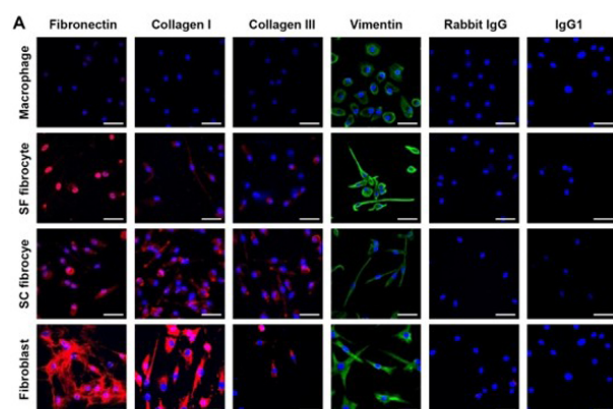


Immunofluorescence Microscopy

Fibrocytes were generated from PBMC under serum-free (SF) and serum-containing (SC) culture conditions, as well as macrophages. Fibroblast cell lines were included as a positive control for stromal cell markers. Immunostaining is shown in red for the stromal markers fibronectin, collagens I and III and their rabbit IgG control, with vimentin and its mouse IgG1 control in green (A). Nuclear staining is shown in blue. Bar represents 50 μ m. Data are representative of at least two separate experiments. Fig 2. PMID: 20305780

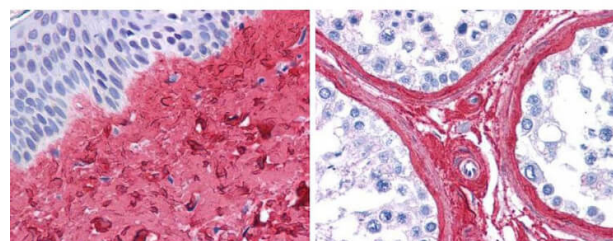


Immunohistochemistry

Rockland anti collagen III antibody (600-401-105 Lot 26016, 1:400, 45 min RT) showed strong staining in FFPE sections of human skin(left, dermis) with moderate to strong red staining and testis (right) where strong staining was observed within connective tissue between seminiferous tubules. The antibody showed strong extracellular staining within connective tissues across many organs with minimal background staining. Slides were steamed in 0.01 M sodium citrate buffer, pH 6.0 at 99-100°C - 20 minutes for antigen retrieval. Images provided courtesy of LifeSpan Biosciences, Seattle, WA

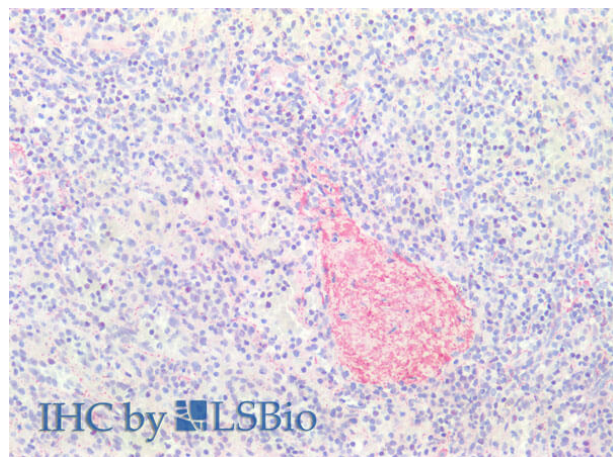


Immunohistochemistry

Immunohistochemistry of Rabbit Anti-Collagen Type III Antibody. Tissue: FFPE normal human spleen tissues (10X). Antigen Retrieval: 0.01 M sodium citrate buffer for 20 minutes. Primary Antibody: Anti-Collagen Type III at 5 μ L/mL for 45 mins at RT. Staining: Anti-Rabbit biotinylated secondary antibody for 30 min at RT. Alkaline phosphatase streptavidin for 30 min at RT. Alkaline phosphatase chromogen substrate for 30 min at RT. The stained slides were evaluated by a pathologist to confirm staining specificity.

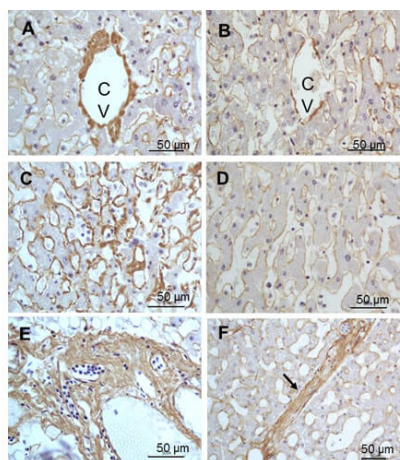

Immunohistochemistry

Immunohistochemistry of Rabbit Anti-collagen type III antibody. Tissue: right lobe of the liver section. A: Central Vein (CV) fibrosis, B: Non-fibrotic CV, C: Perisinusoidal fibrosis, D: Non-fibrotic area, E: Prototypic tract fibrosis, F: Septal fibrosis (arrow). Fixation: formalin fixed paraffin embedded. Antigen retrieval: not required. Primary antibody: Anti-collagen type III at 1:500 for 4°C for 24hr. Secondary antibody: Peroxidase biotin-streptavidin rabbit secondary antibody at 1:10,000 for 45 min at RT. Localization: Anti-collagen type III is intra and extracellular. Staining: 3,3'-diaminobenzidine tetrahydrochloride was used as the chromogen. Nuclei were counterstained purple with hematoxylin.

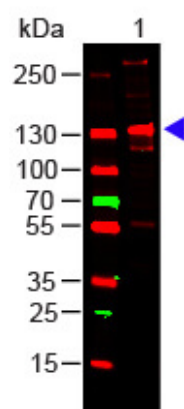

Western Blot

Western Blot of Rabbit Anti-COLLAGEN III Antibody. Lane 1: Human Collagen III (p/n 009-001-105). Load: 100ng per lane. Primary antibody: Collagen III Antibody at 1:1000 o/n at 4°C. Secondary antibody: DyLight™ 649 Goat anti-rabbit (p/n 611-1302) at 1:20,000 for 30 min at RT. Block: MB-070 for 30 min at RT. Predicted/Observed size: 138 kDa, 138 kDa.

References

- Mills SJ et al. Beta-Glucan: A Macrophage Differentiation Modulator That Accelerates the Wound Healing Response. *Int Wound J.* (2026)
- Poret JM et al. Chronic binge alcohol dysregulates omental adipose tissue extracellular matrix in simian immunodeficiency virus-infected macaques. *Alcohol Clin Exp Res (Hoboken).* (2025)
- Baldavira CM et al. Matrisome analysis of NSCLC unveils clinically-important cancer-associated extracellular matrix changes. *Biochim Biophys Acta Mol Basis Dis.* (2025)
- Queiroz ZAJ et al. Type V Collagen as a Critical Regulator of Fibrillar Matrix Remodeling in a Murine Model of Systemic Sclerosis. *Cells.* (2025)
- Zhang Z et al. Management of autofluorescence in formaldehyde-fixed myocardium: choosing the right treatment. *Eur J Histochem.* (2023)
- Szekely T et al. Myelofibrosis progression grading based on type I and type III collagen and fibrillin 1 expression boosted by whole slide image analysis. *Histopathology.* (2023)
- Gerckens M et al. Phenotypic drug screening in a human fibrosis model identified a novel class of antifibrotic therapeutics. *Sci Adv.* (2021)
- Dantas WS, Roschel H, Murai IH, et al. Exercise-Induced Increases in Insulin Sensitivity After Bariatric Surgery Are Mediated By Muscle Extracellular Matrix Remodeling. *Diabetes.* (2020)
- Thomas HM, Ahangar P, Hofma BR, et al. Attenuation of Flightless I Increases Human Pericyte Proliferation, Migration and Angiogenic Functions and Improves Healing in Murine Diabetic Wounds. *Int J Mol Sci.* (2020)
- Tiwari A et al. Loss of HIF1A from pancreatic cancer cells increases expression of ppp1r1b and degradation of p53 to promote invasion and metastasis. *Gastroenterology* (2020)
- Keranov S et al. CILP1 as a biomarker for right ventricular maladaptation in pulmonary hypertension. *Eur Respir J.* (2020)
- Struswick XL et al. Systemic Delivery of Anti-integrin α L Antibodies Reduces Early Macrophage Recruitment, Inflammation, and Scar Formation in Murine Burn Wounds. *Adv Wound Care (New Rochelle).* (2020)
- Kishimoto Y, Yamashita M, Wei A, et al. Reversal of Vocal Fold Mucosal Fibrosis Using siRNA against the Collagen-Specific Chaperone Serpinh1. *Mol Ther Nucleic Acids.* (2019)
- Hikage F, Atkins S, Kahana A, Smith TJ, Chun TH. HIF2A-LOX Pathway Promotes Fibrotic Tissue Remodeling in Thyroid-Associated Orbitopathy. *Endocrinology.* (2019)
- Schulz MC et al. Epithelial-fibroblast cross talk aggravates the impact of the nephrotoxin ochratoxin A. *Biochem Biophys Acta Mol Cell Res.* (2019)
- Noskovičová N, Heinzelmann K, Burgstaller G, Behr J, Eickelberg O. Cub domain-containing protein 1 negatively regulates TGF- β signaling and myofibroblast differentiation. *Am J Physiol Lung Cell Mol Physiol.* (2018)
- Laternser S, Keller H, Leupin O, Rausch M, Graf-Hausner U, Rimann M. A Novel Microplate 3D Bioprinting Platform for the Engineering of Muscle and Tendon Tissues. *SLAS Technol.* (2018)
- Morales MG et al. The pro-fibrotic connective tissue growth factor (CTGF/CCN2) correlates with the number of necrotic-regenerative foci in dystrophic muscle. *J Cell Commun Signal.* (2018)
- Riquelme-Guzmán C, Contreras O, Brandan E. Expression of CTGF/CCN2 in response to LPA is stimulated by fibrotic extracellular matrix via the integrin/FAK axis. *Am J Physiol Cell Physiol.* (2018)

- Li J, Yousefi K, Ding W, Singh J, Shehadeh LA. Osteopontin RNA aptamer can prevent and reverse pressure overload-induced heart failure. *Cardiovasc Res.* (2017)
- Anjum F, Agabalyan NA, Sparks HD, Rosin NL, Kallos MS, Biernaskie J. Biocomposite nanofiber matrices to support ECM remodeling by human dermal progenitors and enhanced wound closure. *Sci Rep.* (2017)
- Rosin NL, Agabalyan N, Olsen K, et al. Collagen structural alterations contribute to stiffening of tissue after split-thickness skin grafting. *Wound Repair Regen.* (2016)
- Budnick I, Hamburg-Shields E, Chen D, et al. Defining the identity of mouse embryonic dermal fibroblasts. *Genesis.* (2016)
- Pirhonen J, Arola J, Sädevirta S, et al. Continuous Grading of Early Fibrosis in NAFLD Using Label-Free Imaging: A Proof-of-Concept Study. *PLoS One.* (2016)
- [View More ...](#)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.