

Datasheet for 600-401-252S**hTERT Antibody****Overview**

| | |
|----------------------|--|
| Description: | Anti-Telomerase catalytic subunit (RABBIT) Antibody - 600-401-252S |
| Item No.: | 600-401-252S |
| Size: | 25 µL |
| Applications: | ELISA, IF, IP, WB, CHIP, FISH, IHC, Multiplex |
| Reactivity: | Human |
| Host Species: | Rabbit |

Product Details

| | |
|----------------------|--|
| Background: | Telomerase is a reverse transcriptase that adds telomeric repeats (TTAGGG) _n to chromosomal ends, compensating for the telomere shortening that occurs with DNA replication. In normal human somatic cells, telomerase is repressed and telomeres progressively shorten, leading to limited lifespan and senescence. Reactivation of telomerase activity is associated with human cancer and cell immortalization. Approximately 85% of human cancers, including breast, prostate, stomach, bladder, colon, and liver cancer, have telomerase activity, whereas most normal somatic cells do not. The specificity of telomerase to human cancer has led to investigations of telomerase activity and expression as a tumor marker. For example, the presence of telomerase activity in human urine has been identified as a marker for human bladder carcinoma. Human telomerase consists of three major subunits: a catalytic protein subunit called hTERT (for human Telomerase Reverse Transcriptase), a template RNA called hTR, and telomerase-associated protein (TEP-1). TERT and hTR are minimally required to reconstitute telomerase activity in vitro. In human cells, hTR is constitutively expressed. TERT transcription is a primary mechanism for regulation of telomerase activity. |
| Synonyms: | rabbit anti-TERT antibody, rabbit anti-Telomerase catalytic subunit antibody, hTERT, Telomerase reverse transcriptase, HEST2, Telomerase-associated protein 2, TP2, EST2, TCS1, TRT |
| Host Species: | Rabbit |
| Clonality: | Polyclonal |
| Format: | IgG |

Target Details

| | |
|-------------------|------|
| Gene Name: | TERT |
|-------------------|------|

| | |
|----------------------------|---|
| Reactivity: | Human |
| Immunogen Type: | Conjugated Peptide |
| Immunogen: | This affinity purified antibody was prepared from whole rabbit serum produced by repeated immunizations with a synthetic peptide corresponding to a region near the carboxy terminal end of hTERT (accession number AF018167). |
| Purity/Specificity: | Affinity Purified Anti-hTERT Antibody was prepared from monospecific antiserum by immunoaffinity chromatography using synthetic peptide coupled to agarose beads. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Rabbit Serum. Although it has been reported that this antibody reacts with mouse TERT (mTERT) (see Drissi, et al. 2001), the binding to mTERT is considerably weaker and less specific than the binding to hTERT (not shown). |
| Relevant Links: | <ul style="list-style-type: none">• 600-401-252 SDS• UniProtKB - O14746• NCBI - NP_001180305.1• GeneID - 7015 |

Application Details

| | |
|--------------------------------|--|
| Tested Applications: | ELISA, IF, IP, WB |
| Suggested Applications: | ChIP, FISH, IHC, Multiplex (Based on references) |

Application Note: Anti-Telomerase catalytic subunit antibody has been tested for use in ELISA, immunoblotting, immunoprecipitation, and immunofluorescence microscopy. In these assays, the antibody detects ectopically-expressed hTERT and high levels of endogenous hTERT. A SY5Y cell nuclear extract can be used as a positive control. This antibody primarily detects hTERT, but several non-specific bands appear on immunoblots. In immunofluorescence microscopy assays, staining with anti-TERT-16 was specific to the nuclei of cells with ectopic TERT expression. In immunoblot assays, whole cell or nuclear extracts were loaded at a concentration of 100 µg protein per well. A working dilution of 1:500 anti-TERT antibody was used followed by a 1:3,000 dilution of HRP goat anti-rabbit IgG as the secondary antibody. For immunofluorescence microscopy staining, a working dilution of 1:500 was used followed by a 1:200 dilution of rhodamine-conjugated donkey anti-rabbit IgG as a secondary antibody. Immunoprecipitation was performed using 20 µL of protein A beads and 2 µL of the anti-TERT serum per 1mg protein from cell lysate. A working dilution of 1:500 is also suggested for immunohistochemistry. To detect TERT, fix cells in 2% paraformaldehyde (in PBS) for 10'. Wash the slides twice in PBS for 5' each. Permeabilize the cells in 0.5% NP-40 for 10'. Wash as before in PBS. Block the cells using PBG buffer (0.2% cold water fish gelatin (Sigma G-7765) and 0.5% BSA in PBS) for 20' at room temperature. Incubate in primary antibody (diluted in PBG) for 1-2 hours at RT or overnight at 4°C. Wash the slides three times in PBG for 5' each. Incubate with secondary antibody (diluted in PBG) for 1 hour at RT in the dark. Wash the slides three times in PBG for 5' each. Mount in DAPI-containing medium.

| | |
|-------------------------|---|
| Assay Dilutions: | All assays should be optimized by the user. Recommended dilutions (if any) may be listed below. |
| ELISA: | 1:10,000 - 1:50,000 |
| IF: | 1:500 |
| IHC: | 1:500 |
| IP: | IP 2µL per mg lysate |
| WB: | 1:500 |

Formulation

| | |
|------------------------|--|
| Physical State: | Liquid (sterile filtered) |
| Concentration: | 1.05 mg/mL by UV absorbance at 280 nm |
| Buffer: | 0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2 |
| Preservative: | 0.01% (w/v) Sodium Azide |
| Stabilizer: | None |


Shipping & Handling

| | |
|----------------------------|---------|
| Shipping Condition: | Dry Ice |
|----------------------------|---------|

Storage Condition: Store vial at -20° C or below prior to opening. This vial contains a relatively low volume of reagent (25 µL). To minimize loss of volume dilute 1:10 by adding 225 µL of the buffer stated above directly to the vial. Recap, mix thoroughly and briefly centrifuge to collect the volume at the bottom of the vial. Use this intermediate dilution when calculating final dilutions as recommended below. Store the vial at -20°C or below after dilution. Avoid cycles of freezing and thawing.


Expiration: Expiration date is three (3) months from date of receipt.

Images

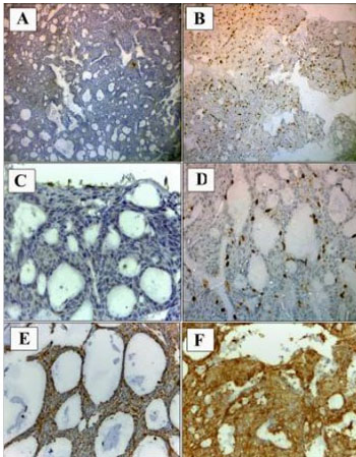


Immunohistochemistry

Immunohistochemistry of proExTMC shows (A) negative expression and (B) positive nuclear expression in HMSC ($\times 100$). High magnification ($\times 400$) images of (C) negative expression and (D) positive nuclear expression of proExTMC in HMSC. Immunohistochemistry of human telomerase reverse transcriptase (hTERT) shows (E) negative expression and (F) positive cytoplasmic expression in HMSC ($\times 400$). Fig 2. PMID: 38584650.

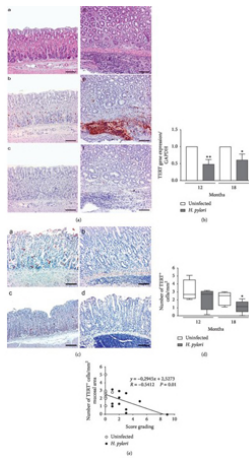


Immunohistochemistry

H. pylori infection decreases mTERT expression in the gastric mucosa of C57BL/6 mice, in the presence of large B lymphocyte aggregates. (a, d) H&E staining and immunostaining of B (b, e) and T (c, f) lymphocytes in gastric sections in infected mice, 12 months after *H. pylori* SS1 infection (d, e, f) and in control mice (a, b, c). Inflammatory infiltrates were observed in the stomach of infected mice, in the lamina propria and submucosa (c). High number of B lymphocytes (e) and a low number of T lymphocytes (f) were present in the inflammatory cell infiltrates in the infected gastric submucosa (e) compared to uninfected (b and c, respectively). Sections of the stomach from the uninfected mice were negative for both B (b) and T (c) lymphocyte staining. Original magnification $\times 4$, bar: 250 μm (a, b, c), and $\times 10$, bar: 100 μm (d, e, f). (b) mTERT gene expression in gastric tissues of *H. pylori* SS1-infected mice at 12 and 18 months after infection quantified by real-time qPCR (Taqman). Results are expressed as means \pm SD of three independent experiments (infected versus uninfected $\square p < 0.05$; $\square\square p < 0.01$). (c) TERT immunolabeling in gastric tissue sections from uninfected mice (a, c) and *H. pylori* SS1-infected (b, d) after 12 (a, b) and 18 (c, d) months. Lower TERT staining is observed in the gastric mucosa in the area of the inflammatory B lymphocyte infiltrates in infected samples. Original magnification: $\times 10$, bar: 100 μm (a, b), and $\times 4$, bar: 250 μm (c, d). (d) Number of TERT-positive cells/mm² mucosal area in gastric tissue sections of uninfected and infected samples at 12 and 18 months. Results are expressed as mean \pm SD (infected versus uninfected $\square p < 0.05$) according to Mann–Whitney analysis. (e) Inverse correlation between the number of TERT-positive cells/mm² mucosal area and the total score grading inflammatory lesions in uninfected (white symbols) and infected mice (black symbol), indicating that TERT level decreases with the exacerbation of gastric inflammation. Each symbol represents one mouse. Fig 3. PMID: 32082377

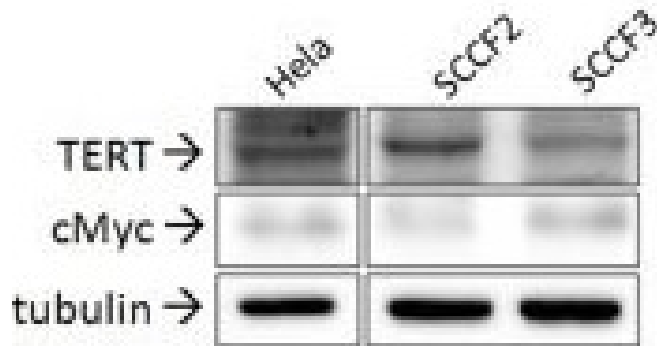
C

Western Blot

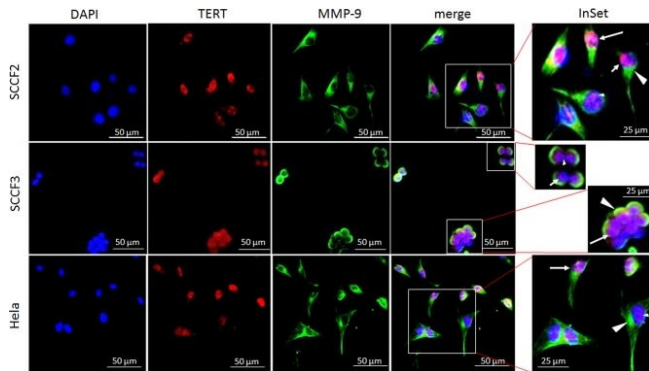Expression of TERT, cMyc, and telomerase activity in SCCF2 and SCCF3 (A,B) TERT and cMyc gene expression in SCCF2 and SCCF3. Data are normalized for β 2-microglobulin expression as housekeeping gene and expressed as relative quantization by the $2^{-\Delta\Delta Ct}$ method. Relative mRNA levels of TERT were lower while cMyc gene expression was higher in SCCF3 with respect to SCCF2 in at least three repeated, independent experiments. (C) A representative Western blotting (WB) gel showing lower TERT protein levels but higher cMyc amounts in SCCF3 vs. SCCF2 is illustrated. HeLa whole cell lysate run along with feline samples confirmed the identity of the band. The blot was stripped and reprobed for tubulin to ensure comparable protein loading and allow normalization. Boxes are cut from the same gel at the same exposure time and properly aligned according to molecular standards loaded onto the gel. Full scans from original gels are shown in Supplementary Figure 2. (D,E) Mean densitometric values \pm SD from at least three repeated, independent WB experiments, normalized for tubulin expression (statistically significant, * $P < 0.05$ and ** $P < 0.01$). (F) Telomerase activity was detected in SCCF2 and SCCF3 by telomeric repeat amplification protocol (TRAP) assay. HeLa cell lysate was run along with feline samples as positive control. A representative gel out of four independent experiments, showing lower telomerase activity in SCCF3 vs. SCCF2 is illustrated (C-: negative control, sample with no lysate; L100bp: 100 base pairs DNA ladder, the first band from the bottom is 100 bp). (G) Quantification of telomerase activity in SCCF2 and SCCF3 by densitometric analysis expressed in % relative to SCCF2. Data were calculated as the ratio between TERT products ladder and 36 bp internal standard and represent the mean \pm standard deviations (SD) of four repeated, independent experiments (statistically significant, ** $P < 0.01$). Figure provided by CiteAb. Source: Front Vet Sci, PMID: 32292795.



Immunocytochemistry

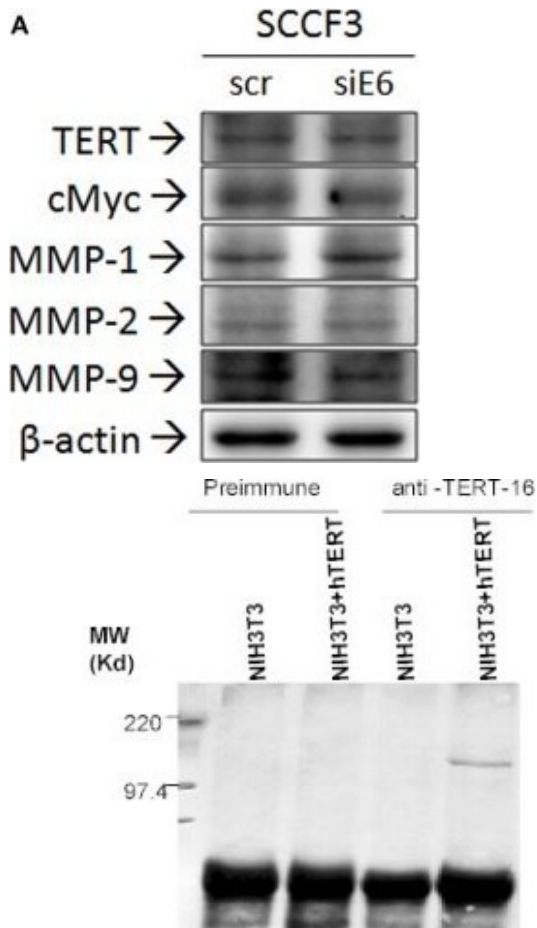Sub-cellular localization of TERT and MMP-9 in SCCF2 and SCCF3. Cells were grown on coverslips and subjected to double indirect IF staining for TERT (red fluorescence) and MMP-9 (green fluorescence). Nuclei were counterstained with DAPI. Inset shows higher magnification of merge panel. TERT was localized in the whole nuclear area (long arrow), compartmentalized in proximity of nuclear membrane (short arrow) or in dot-like spots (small arrowhead). MMP-9 was expressed in the cytoplasm (large arrowhead). HeLa were stained to ensure antibody reactivity. Figure provided by CiteAb. Source: Front Vet Sci, PMID: 32292795.




Western Blot

Effects of FcaPV-2 E6 gene silencing on TERT, cMyc, and MMP expression in SCCF3 cell line. (A) Representative WB panels showing no changes in TERT, cMyc, MMP-1, MMP-2, and MMP-9 protein levels upon FcaPV-2 E6 gene knock-down siRNA (siE6) vs. scramble RNA (scr) treatment. The blot was stripped and reprobred for β-actin to ensure comparable protein loading and allow normalization. (B) Mean densitometric values normalized for β-actin expression ± standard deviations out of two independent experiments. Figure provided by CiteAb. Source: Front Vet Sci, PMID: 32292795.

Western Blot

Immunoprecipitation of hTERT protein from NIH 3T3 cell lysates. The anti-hTERT antibody was used for both immunoprecipitation and western blotting. Anti-hTERT antibody was able to immunoprecipitate TERT protein from cells with ectopic hTERT expression (lane 4). The preimmune serum was unable to immunoprecipitate TERT protein (lanes 1 and 2).

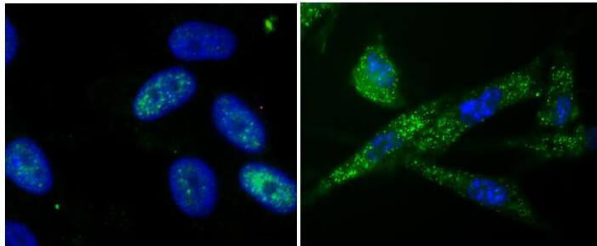


lot 21422, 1:2000 overnight

lot 25694, 1:2000 overnight

Immunofluorescence Microscopy

Immunofluorescence of Rockland anti-hTERT antibody - Two different lots of Rockland anti hTERT antibody were used to stain hTERT on hTERT-over-expressing fibroblasts. Cells were fixed in 4% PFA (in PBS) for 10 min and frozen in -80 after 3 min air-drying before incubation with Rockland anti hTERT 1:2000 overnight and staining with a 1:2000 dilution of Alexafluor488 secondary Ab. Confocal images provided by G. Saretzki, Institute for Ageing and Health, Newcastle University, UK. PMID: 18334557 for more information.

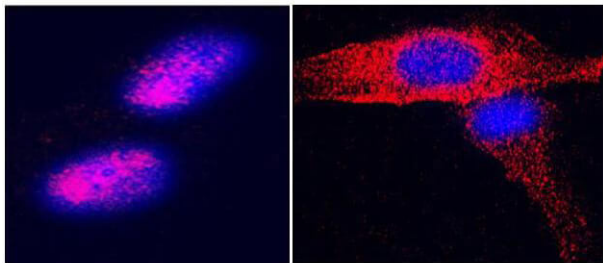


untreated

500 mM H2O2


Immunofluorescence Microscopy

Immunofluorescence of Rockland anti-hTERT antibody was used to stain hTERT on hTERT-over-expressing fibroblasts. Cells were untreated (Left) or treated (Right) with 500 uM H2O2, fixed in 4% PFA (in PBS) for 10 min and frozen in -80 after 3 min air-drying before staining with Rockland anti-hTERT 1:2000 overnight. Confocal images provided by G. Saretzki, Institute for Ageing and Health, Newcastle University, UK. PMID: 18334557 for more information.

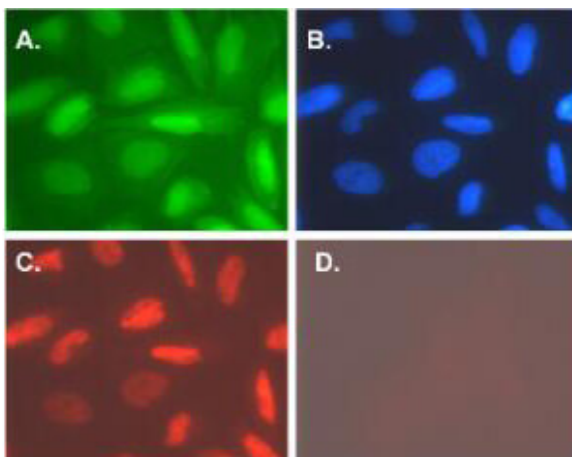


Immunofluorescence Microscopy

Immunofluorescence microscopy of Saos-2 cells transduced with a retroviral vector expressing hTERT and green fluorescent protein (GFP) from an internal ribosomal entry site (IRES). Panel A shows native GFP expression (green), Panel B shows DAPI staining of chromosomes (blue), Panel C shows anti-hTERT staining at a 1:500 dilution followed by washes and addition of a 1:1000 dilution of rhodamine Goat anti-Rabbit IgG (code 611-1002) for detection. Panel D shows no staining of hTERT-transduced cells using pre-immune serum.

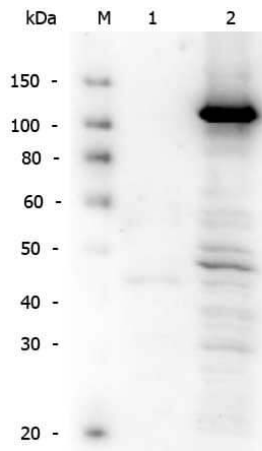
Western Blot

Western Blot of Rabbit anti-Telomerase catalytic subunit antibody. Lane 1: HeLa HV. Lane 2: HeLa 6-2. Load: 10 µg per lane. Primary antibody: Telomerase catalytic subunit antibody at 1:1,000 for overnight at 4°C. Secondary antibody: Peroxidase rabbit secondary antibody (p/n 611-103-122) at 1:40,000 for 45 min at RT. Block: Blocking Buffer for Fluorescent Western Blotting (p/n MB-070). Predicted/Observed size: 127 kDa, 127 kDa for Telomerase catalytic subunit.

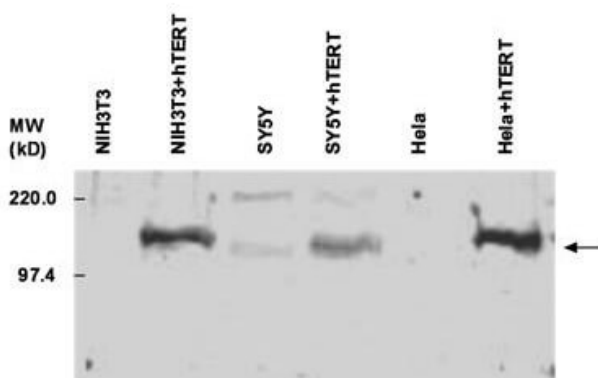
Western Blot

Western blot using anti-hTERT antibody (1:500 dilution). Endogenous levels of mTERT in NIH 3T3 cells (lane 1) and hTERT in HeLa cells (lane 5) are not detectable. After transduction with an hTERT expression vector, both cell types show high levels of hTERT protein (lanes 2 and 6). SY5Y cells, which have high endogenous levels of hTERT, have detectable hTERT protein in both untransduced (lane 3) and transduced (lane 4) cells. The arrow indicates a molecular weight of approximately 127kD, the expected size of hTERT protein.

References

- Brandt A et al. UV1 vaccination in pembrolizumab-treated patients with recurrent or metastatic head and neck cancer: A randomized multicenter phase 2 trial. *Med.* (2025)
- Alabiad MA. et al. Evaluation of Some Prognostic Biomarkers in Human Papillomavirus-Related Multiphenotypic Sinonasal Carcinoma. *Iran J Med Sci.* (2024)
- Muranyi, W et al. Immortalized human choroid plexus endothelial cells enable an advanced endothelial-epithelial two-cell type in vitro model of the choroid plexus. *IScience* (2022)
- Gómez-Ferrer, M et al. HIF-1α and Pro-Inflammatory Signaling Improves the Immunomodulatory Activity of MSC-Derived Extracellular Vesicles. *International Journal of Molecular Sciences* (2021)
- Lagnado, A et al. Neutrophils induce paracrine telomere dysfunction and senescence in ROS-dependent manner. *The Embo Journal* (2021)

- Chatterjee, S et al. Telomerase therapy attenuates cardiotoxic effects of doxorubicin. *Molecular Therapy : the Journal of the American Society of Gene Therapy* (2021)
- Reilly CR et al. The clinical and functional effects of TERT variants in myelodysplastic syndrome. *Blood*. (2021)
- Gissi, C et al. Extracellular vesicles from rat-bone-marrow mesenchymal stromal/stem cells improve tendon repair in rat Achilles tendon injury model in dose-dependent manner: A pilot study. *PLoS One* (2020)
- Altamura G et al. Telomerase Reverse Transcriptase (TERT) Expression, Telomerase Activity, and Expression of Matrix Metalloproteinases (MMP)-1/-2/-9 in Feline Oral Squamous Cell Carcinoma Cell Lines Associated With *Felis catus* Papillomavirus Type-2 Infection. *Front Vet Sci*. (2020)
- Gissi C. et al. Extracellular vesicles from rat-bone-marrow mesenchymal stromal/stem cells improve tendon repair in rat Achilles tendon injury model in dose-dependent manner: A pilot study. *PLoS One* (2020)
- Bussiere FI et al. DNA hypermethylation downregulates telomerase reverse transcriptase (TERT) during *H. pylori*-induced chronic inflammation. *J Oncol*. (2019)
- Perera, ON et al. Telomerase promotes formation of a telomere protective complex in cancer cells. *Science Advances* (2019)
- Delfarah, A et al. Inhibition of nucleotide synthesis promotes replicative senescence of human mammary epithelial cells. *The Journal of Biological Chemistry* (2019)
- Hugdahl et al. Prognostic impact and concordance of TERT promoter mutation and protein expression in matched primary and metastatic cutaneous melanoma. *British Journal of Cancer* (2018)
- Schmidt et al. Dynamics of human telomerase recruitment depend on template- telomere base pairing. *Molecular Biology of the Cell* (2018)
- Kabacik et al. Epigenetic ageing is distinct from senescence-mediated ageing and is not prevented by telomerase expression. *Aging* (2018)
- Bernal et al. Generation of Immortalised But Unstable Cells after hTERT Introduction in Telomere-Compromised and p53-Deficient vHMECs. *International Journal of Molecular Sciences* (2018)
- Li et al. The Werner Syndrome Helicase Coordinates Sequential Strand Displacement and FEN1-Mediated Flap Cleavage during Polymerase δ Elongation. *Molecular and Cellular Biology* (2017)
- Qing H et al. Differential regulation of telomerase reverse transcriptase promoter activation and protein degradation by histone deacetylase inhibition. *J Cell Physiol*. (2016)
- Eitan et al. Excitotoxic and Radiation Stress Increase TERT Levels in the Mitochondria and Cytosol of Cerebellar Purkinje Neurons. *The Cerebellum* (2016)
- Hong et al. Telomerase activates transcription of cyclin D1 gene through an interaction with NOL1. *Journal of Cell Science* (2016)
- Vanderstraeten et al. In Vitro Assessment of the Expression and T Cell Immunogenicity of the Tumor-Associated Antigens BORIS, MUC1, hTERT, MAGE-A3 and Sp17 in Uterine Cancer. *International Journal of Molecular Sciences* (2016)
- Jeong et al. Akt-mediated phosphorylation increases the binding affinity of hTERT for importin α to promote nuclear translocation. *Journal of Cell Science* (2015)

- Theurillat, JPP et al. Prostate cancer. Ubiquitylome analysis identifies dysregulation of effector substrates in SPOP-mutant prostate cancer. *Science (New York, N.Y.)* (2014)
- [View More ...](#)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.