

Datasheet for 600-445-386**VSV-G Epitope Tag Antibody DyLight™ 800 Conjugated****Overview**

Description:	Anti-VSV-G Epitope Tag (RABBIT) Antibody DyLight™ 800 Conjugated - 600-445-386
Item No.:	600-445-386
Size:	100 µg
Applications:	WB
Reactivity:	VSV-G-Tag
Host Species:	Rabbit

Product Details

Background:	In order to improve expression levels, solubility, folding, purification and detection of recombinant proteins, a very common strategy is the fusion of peptides or proteins also known as “tags”, to the target protein. Because these tags are entities with known sequences and well characterized physicochemical properties, they are an essential tool in molecular biology that facilitates expression and purification of recombinant proteins. Because fusion tags constitute themselves antigenic epitopes for which antibodies can be developed they particularly useful for specific detection of the target protein. This Anti-VSV-G Epitope Tag Antibody generated in rabbit is conjugated to DyLight™800.
Synonyms:	Rabbit Anti-VSV-G Epitope Tag DyLight 800™ Conjugated Antibody, Rabbit Anti VSV-G Epitope Tag DyLight 800™ Conjugated Antibody, Rabbit Anti-VSV-G Tag Antibody DyLight 800™ Conjugation
Host Species:	Rabbit
Conjugate:	DyLight™ 800
Clonality:	Polyclonal
Format:	IgG
F/P Ratio:	1.7

Target Details

Reactivity:	VSV-G-Tag
Immunogen Type:	Conjugated Peptide

Immunogen:	This affinity purified antibody was prepared from whole rabbit serum produced by repeated immunizations with a synthetic peptide corresponding aa 501-511 (Y-T-D-I-E-M-N-R-L-G-K) of vesicular stomatitis virus glycoprotein (VSV-G) conjugated to KLH using maleimide.
Purity/Specificity:	This affinity purified antibody is directed against the VSV-G epitope tag and is useful in determining its presence in over expressed proteins in various assays. The antibody recognizes the VSV-G epitope tag (Tyr-Thr-Asp-Ile-Glu-Met-Asn-Arg-Leu-Gly-Lys) fused to either the amino- or carboxy- termini of targeted proteins in transfected or transformed cells.

Application Details

Tested Applications:	WB
Application Note:	Anti-VSV-G DyLight800 Conjugated Antibody has been tested by western blot. This product is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms. The emission spectra for this DyLight™ conjugate match the principle output wavelengths of most common fluorescence instrumentation.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
FLISA:	>1:20,000
IF:	>1:5,000
WB:	1:10,000-1:25,000

Formulation

Physical State:	Lyophilized
Concentration:	1.0 mg/mL by UV absorbance at 280 nm
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2
Preservative:	0.01% (w/v) Sodium Azide
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free
Reconstitution Volume:	100 µL
Reconstitution Buffer:	Restore with deionized water (or equivalent)

Shipping & Handling

Shipping Condition:	Ambient
----------------------------	---------

Storage Condition: Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.


Expiration: Expiration date is one (1) year from date of receipt.

Images

Diagram

Properties of DyLight™ Conjugates.


Emission	Color	DyLight™ Dye	Ex/Em (nm)	ϵ (M ⁻¹ cm ⁻¹)	Similar Dyes
Blue		405	400/420	30,000	Alexa™ 405, Cascade Blue
Green		488	493/518	70,000	Alexa™ 488, Cy2®, FITC
Yellow		549	560/568	150,000	Alexa™ 546, Alexa 555, Cy3®, TRITC
Red		649	646/674	250,000	Alexa™ 647, Cy5®
Near Infrared		680	682/715	140,000	Alexa™ 680, Cy5.5®, IRDye™ 700
Infrared		800	770/794	270,000	IRDye™ 800

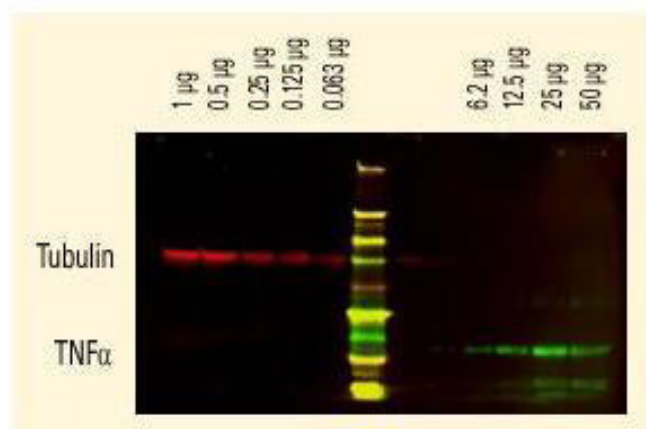

Western Blot

DyLight™ dyes can be used for two-color Western Blot detection with low background and high signal. Anti-tubulin was detected using a DyLight™ 680 conjugate. Anti-TNF α was detected using a DyLight™ 800 conjugate. The image was captured using the Odyssey® Infrared Imaging System developed by LI-COR.

Diagram

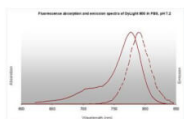

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.