

Datasheet for 611-156-122**Rabbit IgG (H&L) Antibody ATTO 647N Conjugated Pre-Adsorbed****Overview**

Description:	Goat Anti-Rabbit IgG (H&L) Antibody ATTO 647N Conjugated (Min X Bv Ch Gt GP Ham Hs Hu Ms Rt & Sh Serum Proteins) - 611-156-122
Item No.:	611-156-122
Size:	500 µg
Applications:	Dot Blot, WB, IF, IHC, Multiplex
Reactivity:	Rabbit
Host Species:	Goat

Product Details

Background:	Anti-Rabbit IgG (H&L) ATTO 647N Antibody generated in goat detects reactivity to Rabbit IgG. Secreted as part of the adaptive immune response by plasma B cells, immunoglobulin G constitutes 75% of serum immunoglobulins. Immunoglobulin G binds to viruses, bacteria, as well as fungi and facilitates their destruction or neutralization via agglutination (and thereby immobilizing them), activation of the complement cascade, and opsonization for phagocytosis. The whole IgG molecule possesses both the F(c) region, recognized by high-affinity Fc receptor proteins, as well as the F(ab) region possessing the epitope-recognition site. Both the Heavy and Light chains of the antibody molecule are present. Secondary Antibodies are available in a variety of formats and conjugate types. When choosing a secondary antibody product, consideration must be given to species and immunoglobulin specificity, conjugate type, fragment and chain specificity, level of cross-reactivity, and host-species source and fragment composition.
Synonyms:	Goat anti-Rabbit IgG Antibody ATTO647N Conjugation, Goat anti-Rabbit IgG ATTO 647N Conjugated Antibody
Host Species:	Goat
Specificity:	IgG (H&L)
Conjugate:	ATTO 647N
Clonality:	Polyclonal
Format:	IgG
F/P Ratio:	1.72

Target Details

Reactivity:	Rabbit
Immunogen:	Rabbit IgG whole molecule
Purity/Specificity:	Rabbit IgG (H&L) Antibody ATTO 647N was prepared from monospecific antiserum by immunoaffinity chromatography using Rabbit IgG coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Goat Serum, Rabbit IgG and Rabbit Serum. No reaction was observed against Bovine, Chicken, Goat, Guinea Pig, Hamster, Horse, Human, Mouse, Rat and Sheep Serum Proteins. This antibody will react with heavy chains of rabbit IgG and with light chains of most rabbit immunoglobulins.

Application Details

Tested Applications:	Dot Blot, WB
Suggested Applications:	IF, IHC, Multiplex (Based on references)
Application Note:	Anti-Rabbit IgG (H&L) conjugated to ATTO 647N has been tested by dot blot and western blot and is designed for STED microscopy, FRET, immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms. The emission spectra for this ATTO conjugate matches the principle output wavelengths of most common fluorescence instrumentation.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
FLISA:	>1:20,000
IF:	>1:5,000
WB:	>1:10,000

Formulation


Physical State:	Lyophilized
Concentration:	1 mg/mL by UV absorbance at 280 nm
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2
Preservative:	0.01% (w/v) Sodium Azide
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free

Reconstitution Volume:	500 μ L
Reconstitution Buffer:	Restore with deionized water (or equivalent)

Shipping & Handling


Shipping Condition:	Ambient
Storage Condition:	Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.
Expiration:	Expiration date is one (1) year from date of receipt.

Images

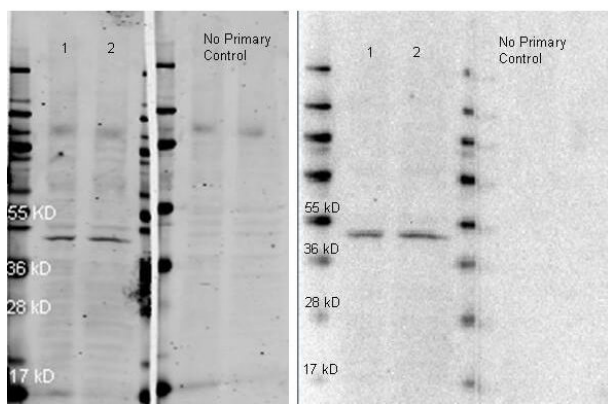


Western Blot

Rockland ATTO 647N conjugated anti rabbit antibody was used to detect anti-Beta Actin antibody (p/n 600-401-886 lot 26928). HeLa (Lane 1) and NIH 3T3 (Lane 2) Whole cell lysates were run on a 4-20% gel, transferred to nitrocellulose under standard conditions, and incubated with anti beta actin at a dilution of 1:2000 (ON 4°C). For secondary antibody detection, blot was incubated for 1 hr RT simultaneously with: 1. Rockland ATTO 647N conjugated anti rabbit antibody (p/n 611-156-122 lot 26426C, 1:10000 in MB-070, Shown on Left and 2. Rockland HRP conjugated anti rabbit IgG (611-1322 lot 19247, 1:10000 in MB-070, shown on right) Blot was dried, imaged at a wavelength of 700 nm on a LiCor Odyssey reader, rewetted in TBS and imaged after 2 min with a 30 sec exposure time using Rockland Femtomax-110 super sensitive Chemiluminescent substrate using the Biorad Versa Doc Imaging System.

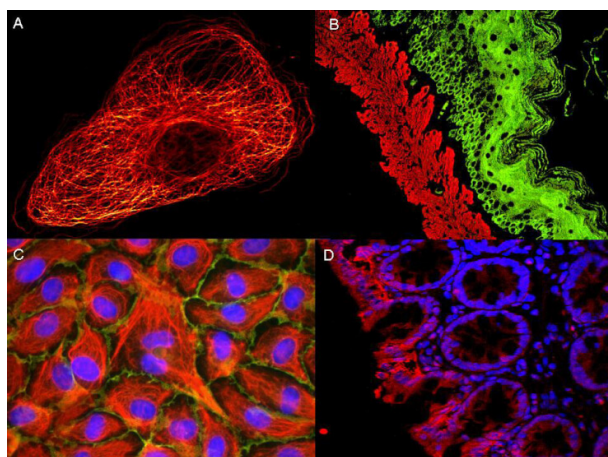

Immunofluorescence Microscopy

ATTO[®] dyes can be used for multicolor immunofluorescent detection with low background and high signal. Examples shown are: A. Tubulin in PtK2- male Rat Kangaroo Kidney Epithelial Cells was detected using ATTO 532 labeled secondary antibody. B. Muscle alpha-actin was stained with a mouse primary antibody and ATTO 488 anti-mouse IgG (green) while Cytokeratin was stained with polyclonal rabbit anti-cytokeratin and ATTO 647N anti-rabbit IgG (red). C. HUVEC (Human umbilical vein endothelial cells) were stained with anti- Vimentin-ATTO 532 (green), anti-E-Cadherin-ATTO 655 (red) and DAPI (blue). D. Rat colon sections were stained with Anti-Aquaporin 3-ATTO 594 antibody. Hoechst 33342 (blue) is used as counterstain. Images provided courtesy of Dr. Jörg Reichwein, ATTO-TEC GmbH

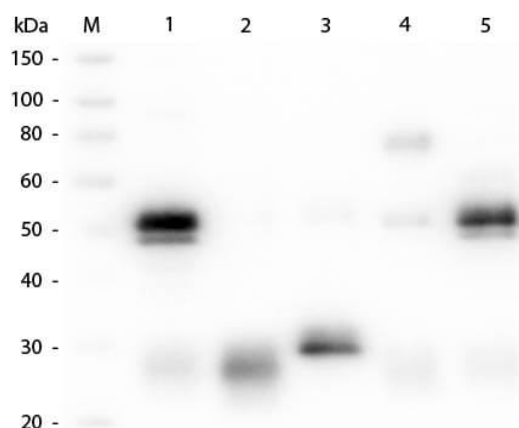

Western Blot

Western Blot of Unconjugated Anti-Rabbit IgG (H&L) (GOAT) Antibody (Min X Bv, Ch, Gt, GP, Ham, Hs, Hu, Ms, Rt & Sh Serum Proteins) (p/n 611-101-122). Lane M: 3 µl Molecular Ladder. Lane 1: Rabbit IgG whole molecule (p/n 011-0102). Lane 2: Rabbit IgG F(ab) Fragment (p/n 011-0105). Lane 3: Rabbit IgG F(c) Fragment (p/n 010-0103). Lane 4: Rabbit IgM Whole Molecule (p/n 011-0107). Lane 5: Normal Rabbit Serum (p/n B309). All samples were reduced. Load: 50 ng per lane. Block: MB-070 for 30 min at RT. Primary Antibody: Anti-Rabbit IgG (H&L) (GOAT) Antibody (Min X Bv, Ch, Gt, GP, Ham, Hs, Hu, Ms, Rt & Sh Serum Proteins) (p/n 611-101-122) 1:1,000 for 60 min at RT. Secondary antibody: Anti-Goat IgG (DONKEY) Peroxidase Conjugated Antibody (p/n CUST10) 1:40,000 in MB-070 for 30 min at RT. Predicted/Observed Size: 25 and 50 kDa for Rabbit IgG and Serum, 25 kDa for F(c) and F(ab), 70 and 23 kDa for IgM. Rabbit F(c) migrates slightly higher.

References

- Shin TW et al. Dense, continuous membrane labeling and expansion microscopy visualization of ultrastructure in tissues. *Nat Commun.* (2025)
- Weissler AL et al. Impaired Presynaptic Function Contributes Significantly to the Pathology of Glycine Receptor Autoantibodies. *Neurol Neuroimmunol Neuroinflamm.* (2025)
- Ogunmowo TH et al. Intersectin and endophilin condensates prime synaptic vesicles for release site replenishment. *Nat Neurosci.* (2025)

- Bao Y et al. A correlative quantitative phase contrast and fluorescence super-resolution microscope for imaging molecules in their cellular context. *bioRxiv [Preprint]*. (2025)
- Coulter O. et al. Adolescent Alcohol Exposure Disrupts Astrocyte-Synaptic Structural And Functional Coupling In The Male Dorsal Hippocampus. *BioRxiv [Preprint]*. (2025)
- Lavogina D et al. Inhibition of epigenetic and cell cycle-related targets in glioblastoma cell lines reveals that onametostat reduces proliferation and viability in both normoxic and hypoxic conditions. *Sci Rep.* (2024)
- Mariani A et al. Neonicotinoid Pesticides Affect Developing Neurons in Experimental Mouse Models and in Human Induced Pluripotent Stem Cell (iPSC)-Derived Neural Cultures and Organoids. *Cells.* (2024)
- Lewis AH et al. Piezo1 ion channels are capable of conformational signaling. *Neuron.* (2024)
- Krause GJ et al. Molecular determinants of the crosstalk between endosomal microautophagy and chaperone-mediated autophagy. *Cell Rep.* (2023)
- Chen J et al. Integrative transcriptomics and cell systems analyses reveal protective pathways controlled by Igfbp-3 in anthracycline-induced cardiotoxicity. *FASEB J.* (2023)
- Pearson-Smith JN et al. Neuronal SIRT3 Deletion Predisposes to Female-Specific Alterations in Cellular Metabolism, Memory, and Network Excitability. *J Neurosci.* (2023)
- Puchov D et al. Analyzing the cellular plasma membrane by fast and efficient correlative STED and platinum replica EM. *Front Cell Dev Biol.* (2023)
- Matthaeus C et al. The molecular organization of differentially curved caveolae indicates bendable structural units at the plasma membrane. *Nat Commun.* (2022)
- Smith, AJ et al. GATD3A, a mitochondrial deglycase with evolutionary origins from gammaproteobacteria, restricts the formation of advanced glycation end products. *Bmc Biology* (2022)
- Hruska M et al. Nanoscale rules governing the organization of glutamate receptors in spine synapses are subunit specific. *Nat Commun.* (2022)
- Cruz-Zaragoza LD et al. An in vitro system to silence mitochondrial gene expression. *Cell.* (2021)
- Scott, WA et al. ATRX proximal protein associations boast roles beyond histone deposition. *PloS Genetics* (2021)
- Kunz TC. et al. Using Expansion Microscopy to visualize and characterize the morphology of mitochondrial cristae. *Front Cell Dev Biol.* (2020)
- Ivanova D et al. CtBP1-mediated membrane fission contributes to effective recycling of synaptic vesicles. *Cell Rep.* (2020)
- Purkey AM et al. AKAP150 palmitoylation regulates synaptic incorporation of Ca²⁺-permeable AMPA receptors to control LTP. *Cell Rep.* (2018)
- Richter KN et al. Comparative synaptosome imaging: a semi-quantitative method to obtain copy numbers for synaptic and neuronal proteins. *Sci Rep.* (2018)
- Hruska et al. Synaptic nanomodules underlie the organization and plasticity of spine synapses. *Nature Neuroscience* (2018)
- Truckenbrodt et al. Newly produced synaptic vesicle proteins are preferentially used in synaptic transmission. *The EMBO Journal* (2018)

- Gomes de Castro et al. Aptamers provide superior stainings of cellular receptors studied under super-resolution microscopy. *PLOS One* (2017)
- Scholefield, J et al. Super-resolution microscopy reveals a preformed NEMO lattice structure that is collapsed in incontinentia pigmenti. *Nature Communications* (2016)
- Saka, SK et al. Correlated optical and isotopic nanoscopy. *Nature Communications* (2014)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.