

Datasheet for B501-0500**Blotto Immunoanalytical Grade (Non-Fat Dry Milk)****Overview**

Description:	Blotto Immunoanalytical Grade (Non-Fat Dry Milk) - B501-0500
Item No.:	B501-0500
Size:	500 g
Applications:	WB, Biochemical Assay, ELISA, IP

Product Details

Background:	Rockland Immunochemicals produces the following blocking buffer reagents: BSA, BLOTTO, BLOTTO A, BLOTTO B and Blocking Buffer for Fluorescent Western Blotting (see related products below). These products are intended to block non-specific binding of proteins in various immunological assays. Typically, western blotting and ELISA methods call for the addition of a blocking agent prior to the addition of primary antibody.
Synonyms:	Blotto Immunoanalytical Grade Non-Fat Dry Milk, Blocking reagent for western, antibody dilution buffer

Target Details

Purity/Specificity:	This product contains highly purified immunoanalytical grade non-fat dry milk.
Relevant Links:	B501-0500 SDS

Application Details

Tested Applications:	WB
Suggested Applications:	Biochemical Assay, ELISA, IP (Based on references)

Application Note: Use Rockland's BLOTTO whenever procedures call for immunoanalytical grade non-fat dry milk. Alternative blocking agents such as BLOTTO A (solution of TBS, 5% non-fat dry milk and 0.05% Tween-20) and BLOTTO B (solution of TBS, 1% non-fat dry milk, 1% BSA and 0.05% Tween-20) are also available. BLOTTO A is a convenient general purpose blocking agent. BLOTTO B is recommended for blocking when phospho-specific antibodies are used. BLOTTO A and BLOTTO B provide Tween-20 in a separate vial in liquid form (p/n TW0025). Users should consider the addition of protease inhibitors and phosphatase inhibitors when required. Typically 50 mM NaF can be added to BLOTTO solutions to inhibit phosphatase activity.

Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
ELISA:	User Optimized
IF:	User Optimized
WB:	User Optimized

Formulation

Physical State:	Lyophilized
Concentration:	1X
Buffer:	None
Preservative:	None
Stabilizer:	None
Reconstitution Buffer:	Restore with deionized water (or equivalent)


Shipping & Handling

Shipping Condition:	Ambient
Storage Condition:	Store container at room temperature prior to opening. After reconstitution, use blocking buffers immediately. Dilute solutions may be stored at 4° C for up to four (4) days. Solutions containing BLOTTO may be frozen.
Expiration:	Expiration date is one (1) year from date of receipt.

Images


Western Blot

Screening multiple cell lines reveals that T47D and MCF7 cells efficiently process recombinant MPO into the mature heterotetrameric form. (A) Cell extracts were made from the

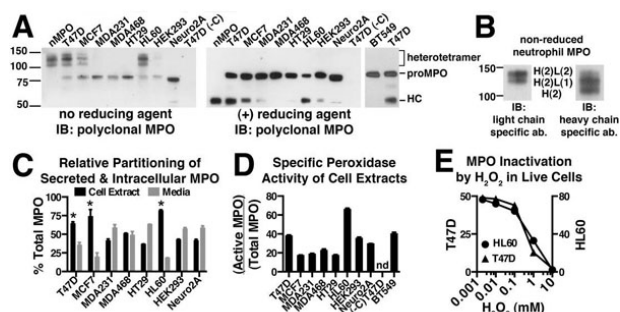

indicated stable cell lines and the total MPO concentration in each extract quantitated by ELISA. Cell extract containing 2 ng of total MPO was run on Bis-Tris gels in the absence or presence of reducing agent and immunoblotted with a multi-epitope polyclonal MPO antibody to assess the MPO species present in each cell line. Under non-reducing conditions the MPO heterotetramer migrates as three bands between the 100–150 kDa markers. Under reducing conditions, the heterotetramer falls apart and we instead observe the heavy chain (HC) of MPO that migrates just above the 50 kDa marker. Monomeric proMPO migrates just above the 75 kDa marker and shifts slightly upward upon reduction of the intra-molecular disulfide bonds.

Untransfected T47D cells (-C) were included as a negative control. (B) The three bands between the 100–150 kDa markers that represent fully processed MPO on non-reducing SDS PAGE were further analyzed to determine if subunit structure was responsible for their distinct gel mobilities. Purified nMPO was run on non-reducing SDS-PAGE and duplicate blots were probed with antibodies specific for either the light chain (L) or heavy chain (H) of MPO. The subunit stoichiometry consistent with the relative immunoreactivity of each band is indicated between the two panels. (C) MPO concentrations in 24 h conditioned media and corresponding whole cell extracts were determined by ELISA and the relative amounts of secreted versus cellular MPO were calculated for each cell line. (D) Whole cell extracts from each cell line were assayed for peroxidase activity and total MPO content and the relative specific activity (Active MPO/Total MPO). (E) Comparison of the sensitivity of native (HL60) and recombinant (T47D) MPO to inactivation by H₂O₂. Live cells were treated with increasing concentrations of H₂O₂ in the growth media for 1.5 h at which time further inactivation was stopped by the addition of catalase. Specific peroxidase activity was measured in cell extracts generated from treated cells as in panel D. IC₅₀ values are 2.8 mM (HL60) and 0.96 mM (T47D-MPO). (C, D & E) Assay points are triplicate measurements and plotted as mean ± SE. The data shown is representative of at least three independent experiments. Membranes were blocked with 4% non-fat dry milk and plates were blocked 2% non-fat dry milk (p/n B501-0500). (nd) None detected. Fig 2. PMID: 26890638.

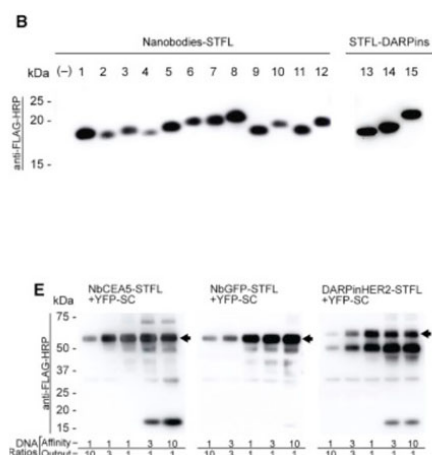


Western Blot

FD-CF Production of On-Demand Affinity Products

(B) Anti-FLAG WB showing FD-CF expression of 12 Nanobodies and three DARPins containing ST-FLAG tags (STFL). Specific antigen targets are: (1) CEA5, (2) dengue Virus NS1, (3) GFP, (4) HIV p23-gag, (5) norovirus capsid VP1, (6) rotavirus capsid VP6, (7) C. difficile exotoxin TcdA, (8) P. falciparum VAR2CSA, (9) GLUT1, (10) mCherry, (11) Vimentin, (12) Glycophorin A, (13) HER-2, (14) VEGF-A, and (15) epCAM. (–) indicates a DNA null control reaction.

(E) Anti-FLAG WB showing one-pot manufacturing of affinity-output proteins, produced by mixing DNA templates encoding different affinity components (left, anti-CEA5-STFL Nanobody; middle, anti-GFP-STFL Nanobody; right, STFL-anti-HER2 DARPin) with DNA encoding the YFP-SC output component in a single FD-CF reaction. Template ratios are shown below. Membrane was blocked in 4% Blotto (p/n B501-0500) + 2% cold water fish gelatin. Figure 4. PMID: 27662092.

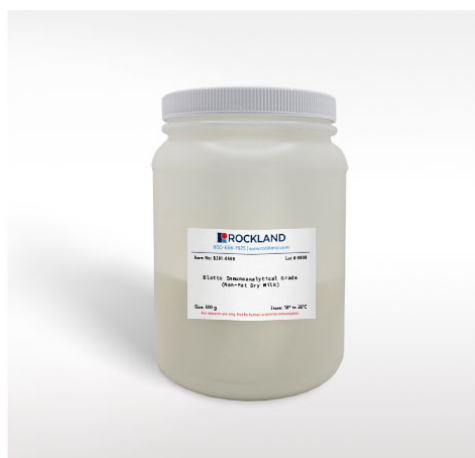

Bottle

Blotto Immunoanalytical Grade (Non-Fat Dry Milk)

References

- Lin M et al. Tetramethylthiuram disulfide induces stress granules and DNA damage through oxidative stress in human lung epithelial cells. *Anim Cells Syst (Seoul)*. (2026)
- Jalkanen J et al. Cytoarchitectural multi-depot profiling reveals immune-metabolic crosstalk in human colon-associated adipose tissue. *Cell Metab*. (2026)
- Takao S et al. Epigenetic mechanisms controlling human leukemia stem cells and therapy resistance. *Nat Commun*. (2025)

- Hahn PA et al. Exogenous estrogen enhances T cell activation in male primates. *Cell Rep.* (2025)
- Kim SE et al. Biological Activities of Crude Grape Tendril Extract Against Oxidative Stress, Inflammation, and Allergy. *Sage.* (2025)
- Chen J et al. Specific macrophage RhoA targeting CRISPR-Cas9 for mitigating osteoclastogenesis-induced joint damage in inflammatory arthritis. *Cell Rep Med.* (2025)
- Joshi R et al. A Novel Bispecific Integrin $\alpha 5\beta 1/\alpha v$ Antibody Reprograms the Myc-Regulated Basal Phenotype of Prostate Cancer with NK Cell-Mediated Tumor Elimination. *Mol Cancer Res.* (2025)
- Hahn PA et al. The Immunogenicity of AAV-Encoded HIV-1 bNAbs in Rhesus Macaques Is Unaffected by a Short Course of the Immunomodulator CTLA4Ig. *AIDS Res Hum Retroviruses.* (2025)
- Iksat N et al. Heat Stress Induces Partial Resistance to Tomato Bushy Stunt Virus in *Nicotiana benthamiana* Via Combined Stress Pathways. *Viruses.* (2025)
- Lee HJ et al. Neuroprotection against beta-amyloid toxicity by the novel estrogen receptor modulator STX requires convergent signaling pathways. *Front Mol Neurosci.* (2025)
- Cho N et al. Alternative splicing of KRAS exon 4 promotes tumor progression via enhanced KRAS4A oncogenic activity. *Anim Cells Syst (Seoul).* (2025)
- Pfeifle A et al. *Borrelia burgdorferi* Strain-Specific Differences in Mouse Infectivity and Pathology. *Pathogens.* (2025)
- Ardeshtir A et al. Determinants of successful AAV-vectored delivery of HIV-1 bNAbs in early life. *Nature.* (2025)
- Kutzler MA et al. The molecular immune modulator adenosine deaminase-1 enhances HIV specific humoral and cellular responses to a native-like HIV envelope trimer DNA vaccine. *Res Sq [preprint]* (2024)
- Euliano EM et al. Intra-lymph node crosslinking of antigen-bearing polymers enhances humoral immunity and dendritic cell activation. *Bioeng Transl Med.* (2024)
- Euliano EM et al. A TLR7 Agonist Conjugated to a Nanofibrous Peptide Hydrogel as a Potent Vaccine Adjuvant. *bioRxiv [preprint]* (2024)
- Moya-Alvarado G et al. PLC- γ -Ca²⁺ pathway regulates axonal TrkB endocytosis and is required for long-distance propagation of BDNF signaling. *Front Mol Neurosci.* (2024)
- Askari F et al. Phosphatidylinositol 3-phosphate regulates iron transport via PI3P-binding CgPIL1 protein. *Cell Rep.* (2023)
- Choi S et al. Polyhexamethylene guanidine phosphate increases stress granule formation in human 3D lung organoids under respiratory syncytial virus infection. *Ecotoxicol Environ Saf.* (2022)
- Ashimova A et al. Sustained Delivery of a Monoclonal Antibody against SARS-CoV-2 by Microencapsulated Cells: A Proof-of-Concept Study. *Pharmaceutics.* (2022)
- Quinn JF et al. The novel estrogen receptor modulator STX attenuates Amyloid- β neurotoxicity in the 5XFAD mouse model of Alzheimer's disease. *Neurobiol Dis.* (2022)
- Aronson, M.R., Dahl, E.S., Halle, J.A. et al. Re-engineering Antimicrobial Peptides into Oncolytics Targeting Drug-Resistant Ovarian Cancers. *Cel. Mol. Bioeng.* (2020)
- Maki RA, Holzer M, Motamedchaboki K, et al. Human myeloperoxidase (hMPO) is expressed in neurons in the substantia nigra in Parkinson's disease and in the hMPO- α -synuclein-A53T mouse model, correlating with increased nitration and aggregation of α -synuclein and exacerbation of motor impairment. *Free Radic Biol Med.* (2019)

- Ferko MA et al. Effects of metal ions on caspase-1 activation and interleukin-1 β release in murine bone marrow-derived macrophages. *PLoS One*. (2018)
- Laura RP et al. T47D Cells Expressing Myeloperoxidase Are Able to Process, Traffic and Store the Mature Protein in Lysosomes: Studies in T47D Cells Reveal a Role for Cys319 in MPO Biosynthesis that Precedes Its Known Role in Inter-Molecular Disulfide Bond Formation. *PLoS One*. (2016)
- Pardee K et al. Portable, on-demand biomolecular manufacturing. *Cell*. (2016)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.