

Datasheet for KOA0302**Mouse P-Selectin ELISA Kit****Overview**

Description:	Mouse P-Selectin AccuSignal ELISA Kit - KOA0302
Item No.:	KOA0302
Size:	1 Kit
Applications:	ELISA
Reactivity:	Mouse
Range:	62.5 pg/ml - 4000 pg/ml

Product Details

Background:	P-selectin, also called GMP-140, CD62, or selectin P, is a 140-kD adhesion molecule, expressed at the surface of activated cells, that mediates the interaction of activated endothelial cells or platelets with leukocytes. It is stored in secretory granules and expressed at the plasma membrane after cell activation. It is known to play an important role in atherosclerosis. The major ligand for P-selectin on leukocytes is P-selectin glycoprotein ligand-1(PSGL-1). The standard product used in this kit is recombinant mouse P-Selectin, excluding intercellular P-Selectin and transmembrane domain. It has 42-709 amino acids sequence with the molecular mass of 99.1KDa. As a result of glycosylation, the molecular mass of 190-200KDa is revealed by SDS-PAGE.
Synonyms:	Antigen CD62, CD62, CD62 antigen like family member P, GMP140, Granule membrane protein 140, Granule membrane protein 140kDa, GRMP, LECAM3, Leukocyte endothelial cell adhesion molecule 3, Leukocyte-endothelial cell adhesion molecule 3, LYAM3_HUMAN, P selectin, P-selectin, PADGEM, Platelet activation dependent granule external membrane protein, Platelet activation dependent granule-external membrane protein, Platelet alpha granule membrane protein, PSEL, Selectin P
Detection Kit Type:	ELISA Kit
Detection Range:	62.5 pg/ml - 4000 pg/ml
Sensitivity:	<5 pg/ml

Target Details

Gene Name:	SELP
-------------------	------

Reactivity:	Mouse
Immunogen:	Expression system for standard: NSO; Immunogen sequence: W42-A709
Purity/Specificity:	Natural and recombinant mouse sP-Selectin. There is no detectable cross-reactivity with other relevant proteins.
Relevant Links:	• UniProtKB - Q01102 • NCBI - NP_035477.1 • GeneID - 20344 • KOA0302 Protocol

Application Details

Tested Applications:	ELISA
Application Note:	Useful in Sandwich ELISA for Quantitative Detection of Antigen. Aliquot 0.1ml per well of the 4000pg/ml, 2000pg/ml, 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml mouse P-Selectin standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Zero well). Add 0.1ml of each properly diluted sample of mouse cell culture supernates, serum or plasma(citrate) to each empty well. It is recommended that each mouse P-Selectin standard solution and each sample be measured in duplicate.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.


Formulation

Concentration:	1 Kit
Anticoagulant:	Sodium Citrate

Shipping & Handling

Shipping Condition:	Wet Ice
Storage Condition:	Store vials at 4°C prior to opening. Centrifuge product if not completely clear after standing at room temperature. This product is stable for 6 months at 4°C as an undiluted liquid. Dilute only prior to immediate use. For extended storage freeze at -20°C or below for 12 months. Avoid cycles of freezing and thawing.
Expiration:	See kit insert for complete instructions.

Images

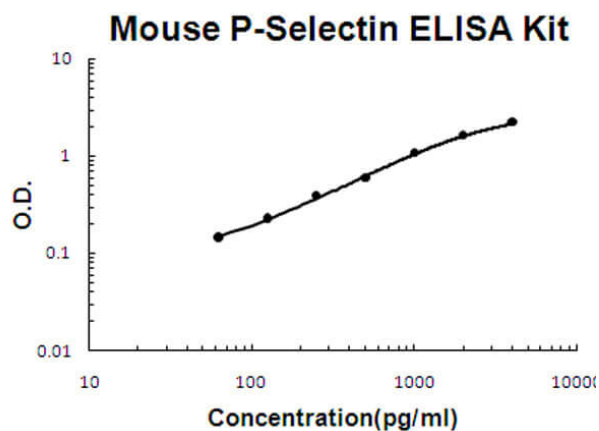

ELISA

Mouse P-Selectin AccuSignal ELISA Kit standard curve. Assay Range: 62.5pg/ml-4000pg/ml. Sensitivity: <5pg/ml.

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.