

Datasheet for KOA0315

Human TIMP-3 ELISA Kit**Overview**

Description:	Human TIMP-3 AccuSignal ELISA Kit - KOA0315
Item No.:	KOA0315
Size:	1 Kit
Applications:	ELISA
Reactivity:	Human
Range:	156 pg/ml - 10000 pg/ml

Product Details

Background:	The tissue inhibitors of metalloproteinases (TIMPs) are natural inhibitors of the matrix metalloproteinases, a group of zinc-binding endopeptidases involved in the degradation of the extracellular matrix. The TIMP3 gene is expressed in many tissues, with highest expression in the placenta. TIMP3 encodes a potent angiogenesis inhibitor and is mutated in Sorsby fundus dystrophy, a macular degenerative disease with submacular choroidal neovascularization. TIMP3 gene is mapped to 22q12.1-q13.2. Mutations in TIMP3 cause the autosomal dominant disorder Sorsby's fundus dystrophy (SFD).
Synonyms:	HSMRK222, K222, K222TA2, Metalloproteinase inhibitor 3, MIG 5 protein, Protein MIG-5, SFD, TIMP metalloproteinase inhibitor 3, TIMP-3, TIMP3, TIMP3_HUMAN, Tissue inhibitor of metalloproteinases 3
Detection Kit Type:	ELISA Kit
Detection Range:	156 pg/ml - 10000 pg/ml
Sensitivity:	<2 pg/ml

Target Details

Gene Name:	TIMP3
Reactivity:	Human
Immunogen:	Expression system for standard: NSO; Immunogen sequence: C24-P211
Purity/Specificity:	Natural and recombinant human TIMP-3. There is no detectable cross-reactivity with other relevant proteins.

Relevant Links:	<ul style="list-style-type: none">• UniProtKB - P35625• NCBI - NP_000353.1• GeneID - 7078• KOA0315 Protocol
------------------------	--

Application Details

Tested Applications:	ELISA
Application Note:	Useful in Sandwich ELISA for Quantitative Detection of Antigen. Aliquot 0.1ml per well of the 10000pg/ml, 5000pg/ml, 2500pg/ml, 1250pg/ml, 625pg/ml, 312pg/ml, 156pg/ml human TIMP-3 standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Zero well). Add 0.1ml of each properly diluted sample of human cell culture supernates, serum, plasma(heparin , EDTA) or saliva to each empty well. It is recommended that each human TIMP-3 standard solution and each sample be measured in duplicate.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
ELISA:	156pg/ml-10000pg/ml
Other:	DOE: 2-8° C - 12/25/2023 -20° C - 06/24/2024


Formulation

Concentration:	1 Kit
Anticoagulant:	Heparin Sodium

Shipping & Handling

Shipping Condition:	Wet Ice
Storage Condition:	Store vials at 4°C prior to opening. Centrifuge product if not completely clear after standing at room temperature. This product is stable for 6 months at 4°C as an undiluted liquid. Dilute only prior to immediate use. For extended storage freeze at -20°C or below for 12 months. Avoid cycles of freezing and thawing.
Expiration:	See kit insert for complete instructions.

Images

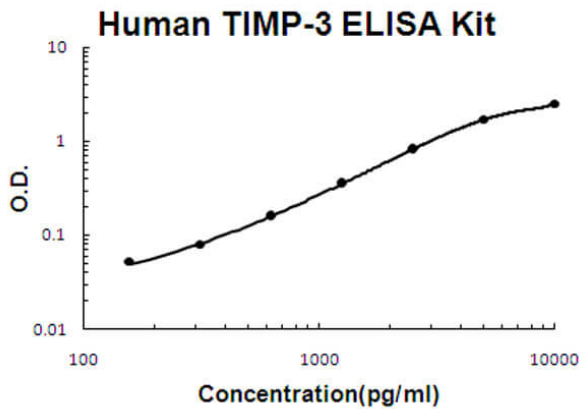

ELISA

Human TIMP-3 AccuSignal ELISA Kit standard curve. Assay Range: 156pg/ml-10000pg/ml. Sensitivity: <2pg/ml.

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.