

Datasheet for KOA0615**Mouse Complement C5a ELISA Kit****Overview**

Description:	Mouse Complement C5a AccuSignal ELISA Kit - KOA0615
Item No.:	KOA0615
Size:	1 Kit
Applications:	ELISA
Reactivity:	Mouse
Range:	15.6 pg/ml - 1000 pg/ml

Product Details

Background:	C5a is a protein fragment released from complement component C5. This gene is mapped to 9q33.2. The protein encoded by this gene is the fifth component of complement, which plays an important role in inflammatory and cell killing processes. This protein is comprised of alpha and beta polypeptide chains that are linked by a disulfide bridge. An activation peptide, C5a, which is an anaphylatoxin that possesses potent spasmogenic and chemotactic activity, is derived from the alpha polypeptide via cleavage with a convertase. The C5b macromolecular cleavage product can form a complex with the C6 complement component, and this complex is the basis for formation of the membrane attack complex, which includes additional complement components. Mutations in this gene cause complement component 5 deficiency, a disease where patients show a propensity for severe recurrent infections. Defects in this gene have also been linked to susceptibility to liver fibrosis and to rheumatoid arthritis.
Synonyms:	anaphylatoxin C5a analog, C3 and PZP-like alpha-2-macroglobulin domain-containing protein 4, C5, CO5_HUMAN, complement C5, Complement C5 alpha chain, Complement component 5, CPAMD4, ECLZB
Detection Kit Type:	ELISA Kit
Detection Range:	15.6 pg/ml - 1000 pg/ml
Sensitivity:	<10 pg/ml

Target Details

Gene Name:	C5
Reactivity:	Mouse

Immunogen:	Expression system for standard: E.coli; Immunogen sequence: N679-R755
Purity/Specificity:	Natural and recombinant mouse Complement C5a. There is no detectable cross-reactivity with other relevant proteins.
Relevant Links:	• SDS • UniProtKB - P06684 • NCBI - NP_034536.2 • GeneID - 15139 • KOA0615 Protocol

Application Details

Tested Applications:	ELISA
Application Note:	Useful in Sandwich ELISA for Quantitative Detection of Antigen. Aliquot 0.1ml per well of the 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml, 31.2pg/ml, 15.6pg/ml mouse Complement C5a standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Zero well). Add 0.1ml of each properly diluted sample of mouse cell culture supernates, serum or plasma(heparin, EDTA) to each empty well. It is recommended that each mouse Complement C5a standard solution and each sample be measured in duplicate.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
ELISA:	15.6pg/ml-1000pg/ml
Other:	DOE: 2-8° C - 11/08/2026 -20° C - 05/08/2027

Formulation

Concentration:	1 Kit
Anticoagulant:	Heparin Sodium


Shipping & Handling

Shipping Condition:	Wet Ice
----------------------------	---------

Storage Condition: Store vials at 4°C prior to opening. Centrifuge product if not completely clear after standing at room temperature. This product is stable for 6 months at 4°C as an undiluted liquid. Dilute only prior to immediate use. For extended storage freeze at -20°C or below for 12 months. Avoid cycles of freezing and thawing.

Expiration: See kit insert for complete instructions.

Images

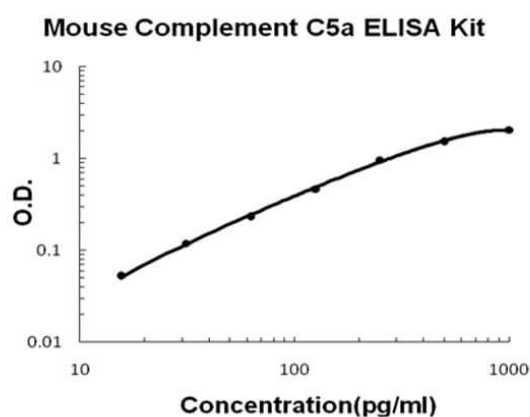

ELISA

Mouse Complement C5a AccuSignal ELISA Kit standard curve. Assay Range: 15.6pg/ml-1000pg/ml. Sensitivity: <10pg/ml.

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.