

Datasheet for KOA0685**Mouse SP-D ELISA Kit****Overview**

Description:	Mouse SP-D AccuSignal ELISA Kit - KOA0685
Item No.:	KOA0685
Size:	1 Kit
Applications:	ELISA
Reactivity:	Mouse
Range:	312 pg/ml - 20,000 pg/ml

Product Details

Background:	Surfactant, pulmonary-associated protein D, also known as SFTPD or SP-D, is a protein which in humans is encoded by the SFTPD gene. By fluorescence in situ hybridization, the SP-D gene was localized in 10q22.2-q23.1. On the basis of homology with other collectins, potential functions for SP-D include roles in innate immunity and surfactant metabolism, SP-D is produced in the bronchiolar and terminal epithelium of human fetal lung from about 21 weeks of gestation. What's more, SP-A and SP-D act as dual-function surveillance molecules that reverse orientation and function and become initiators of host-defense reactions.
Synonyms:	COLEC 7, COLEC7, Collectin-7, Collectin7, Lung surfactant protein D, PSP D, PSP-D, PSP-D Surfactant protein D, PSPD, Pulmonary surfactant apoprotein, Pulmonary surfactant associated protein D, Pulmonary surfactant associated protein PSP-D, Pulmonary surfactant-associated protein D, SFTP 4, SFTP4, SFTPD, SFTPD_HUMAN, SP D, SP-D, Surfactant associated protein pulmonary 4, Surfactant pulmonary associated protein D
Detection Kit Type:	ELISA Kit
Detection Range:	312 pg/ml - 20,000 pg/ml
Sensitivity:	<10 pg/ml

Target Details

Gene Name:	SFTPD
Reactivity:	Mouse
Immunogen:	Expression system for standard: CHO; Immunogen sequence: A20-F374

Purity/Specificity: Natural and recombinant mouse SP-D. There is no detectable cross-reactivity with other relevant proteins.

Relevant Links:

- [UniProtKB - P50404](#)
- [NCBI - NP_033186.1](#)
- [GeneID - 20390](#)
- [KOA0685 Protocol](#)

Application Details

Tested Applications: ELISA

Application Note: Useful in Sandwich ELISA for Quantitative Detection of Antigen. Aliquot 0.1ml per well of the 20,000pg/ml, 10,000pg/ml, 5000pg/ml, 2500pg/ml, 1250pg/ml, 625pg/ml, 312pg/ml mouse SP-D standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Zero well). Add 0.1ml of each properly diluted sample of mouse cell culture supernates, serum or plasma (heparin, EDTA) to each empty well. It is recommended that each mouse SP-D standard solution and each sample be measured in duplicate.

Assay Dilutions: All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.

Other: DOE:
2-8° C - 06/03/2020
-20° C - 12/02/2020

Formulation

Concentration: 1 Kit

Anticoagulant: Heparin Sodium


Shipping & Handling

Shipping Condition: Wet Ice

Storage Condition: Store vials at 4°C prior to opening. Centrifuge product if not completely clear after standing at room temperature. This product is stable for 6 months at 4°C as an undiluted liquid. Dilute only prior to immediate use. For extended storage freeze at -20°C or below for 12 months. Avoid cycles of freezing and thawing.

Expiration: See kit insert for complete instructions.

Images

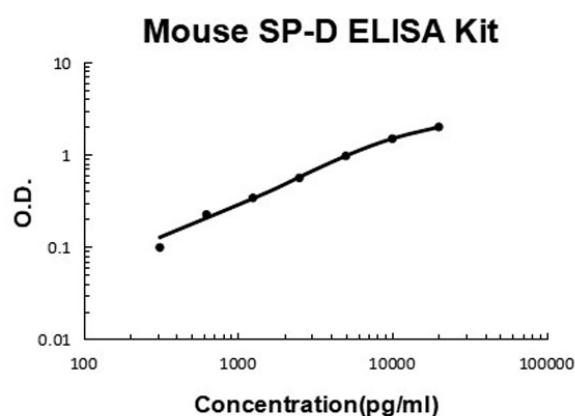

ELISA

Mouse SP-D ELISA Kit standard curve.

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.