

Datasheet for KOA0859**Human CD73 - NT5E ELISA Kit****Overview**

Description:	Human CD73/NT5E AccuSignal ELISA Kit - KOA0859
Item No.:	KOA0859
Size:	1 Kit
Applications:	ELISA
Reactivity:	Human
Range:	31.2 pg/ml - 2000 pg/ml

Product Details

Background:	5'-nucleotidase (5'-NT), also known as ecto-5'-nucleotidase or CD73 (Cluster of Differentiation 73), is an enzyme that in humans is encoded by the NT5E gene. Ecto-5-prime-nucleotidase catalyzes the conversion at neutral pH of purine 5-prime mononucleotides to nucleosides, the preferred substrate being AMP. The enzyme consists of a dimer of 2 identical 70-kD subunits bound by a glycosyl phosphatidyl inositol linkage to the external face of the plasma membrane. And the enzyme is used as a marker of lymphocyte differentiation. Consequently, a deficiency of NT5 occurs in a variety of immunodeficiency diseases. Other forms of 5-prime nucleotidase exist in the cytoplasm and lysosomes and can be distinguished from ecto-NT5 by their substrate affinities, requirement for divalent magnesium ion, activation by ATP, and inhibition by inorganic phosphate. Rare allelic variants are associated with a syndrome of adult-onset calcification of joints and arteries (CALJA) affecting the iliac, femoral, and tibial arteries reducing circulation in the legs and the joints of the hands and feet causing pain.
Synonyms:	5' NT, 5' nucleotidase (CD73) 5' nucleotidase precursor, 5' nucleotidase ecto, 5' nucleotidase ecto (CD73), 5'-NT, 5'-nucleotidase, 5NTD_HUMAN, CD73, CD73 antigen, E5NT, Ecto 5' nucleotidase, Ecto-5'-nucleotidase, eN, eNT, NT, NT5, NT5E, NTE, Purine 5 Prime Nucleotidase
Detection Kit Type:	ELISA Kit
Detection Range:	31.2 pg/ml - 2000 pg/ml
Sensitivity:	<10 pg/ml

Target Details

Gene Name:	NT5E
-------------------	------

Reactivity:	Human
Immunogen:	Expression system for standard: CHO; Immunogen sequence: W27-K547
Purity/Specificity:	Natural and recombinant human CD73. There is no detectable cross-reactivity with other relevant proteins.
Relevant Links:	• UniProtKB - P21589 • NCBI - NP_001191742.1 • GenelD - 4907 • KOA0859 Protocol

Application Details

Tested Applications:	ELISA
Application Note:	Useful in Sandwich ELISA for Quantitative Detection of Antigen. Aliquot 0.1ml per well of the 2000pg/ml, 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml, 31.2pg/ml human CD73 standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Zero well). Add 0.1ml of each properly diluted sample of human cell culture supernates, serum and plasma (heparin, EDTA) to each empty well. It is recommended that each human CD73 standard solution and each sample be measured in duplicate.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
Other:	DOE: 2-8° C - 07/16/2020 -20° C - 01/14/2021

Formulation

Concentration:	1 Kit
Anticoagulant:	Heparin Sodium


Shipping & Handling

Shipping Condition:	Wet Ice
Storage Condition:	Store vials at 4°C prior to opening. Centrifuge product if not completely clear after standing at room temperature. This product is stable for 6 months at 4°C as an undiluted liquid. Dilute only prior to immediate use. For extended storage freeze at -20°C or below for 12 months. Avoid cycles of freezing and thawing.

Expiration: See kit insert for complete instructions.

Images

Human CD73 ELISA Kit

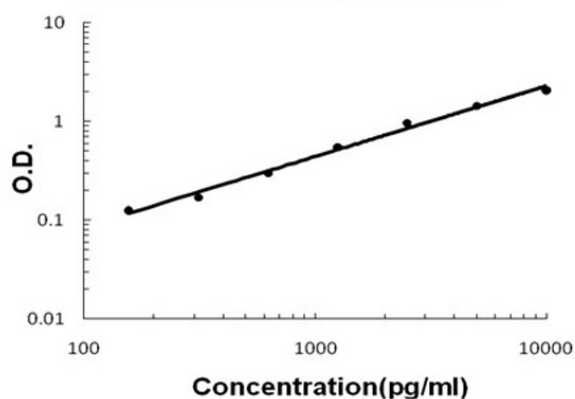

ELISA

Human CD73 ELISA Kit standard curve.

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.