

Datasheet for MB-014**0.5 M EDTA pH 8.0****Overview**

Description:	0.5 M EDTA pH 8.0 - MB-014
Item No.:	MB-014
Size:	100 mL
Applications:	FC

Product Details

Background:	EDTA is widely used for scavenging metal ion in biochemistry and molecular biology. Ion depletion is commonly used to deactivate metal-dependent enzymes, either as an assay for their reactivity or to suppress damage to DNA or proteins. In tissue culture EDTA is used as a chelating agent that binds to calcium and prevents joining of cadherins between cells, preventing clumping of cells grown in liquid suspension, or detaching adherent cells for passaging. EDTA is also known to inhibit a range of metalloproteinases, the method of inhibition occurs via the chelation of the metal ion required for catalytic activity.
Synonyms:	Ethylenediaminetetraacetic acid, EDTA

Target Details

Relevant Links:	MB-014 SDS
------------------------	--

Application Details

Suggested Applications:	FC (Based on references)
Application Note:	Ethylenediaminetetraacetic acid is a concentrated stock solution and should be diluted appropriately with distilled, deionized water or equivalent to its final working concentration. This buffer consists of 0.5 M Trisodium EDTA at pH 8.0. Meticulously prepared using ultra pure reagents dissolved in highly polished pharmaceutical grade deionized water (DI) treated with diethyl pyrocarbonate (DEPC).
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.


Formulation

Physical State:	Liquid (sterile filtered)
Concentration:	0.5M
Buffer:	See application note.
Preservative:	None
Stabilizer:	None

Shipping & Handling

Shipping Condition:	Ambient
Storage Condition:	Store container at room temperature (18° to 26° C) prior to opening. If desired, the solution may be stored at 4° C or less. Some salts may precipitate out of solution at lower temperature. Allow buffer to equilibrate to room temperature (18° to 26° C) to restore solubility of some salts.
Expiration:	Expiration date is six (6) months from date of receipt.


Images

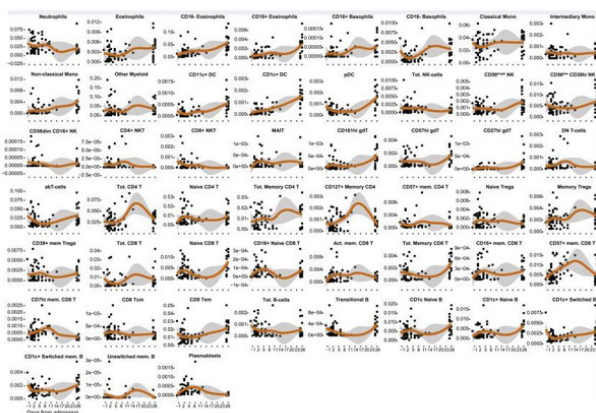

Flow Cytometry

Immune Cell Proportions in COVID-19.

Proportion of 57 white blood cell populations determined by mass cytometry from acute to recovery phase of COVID-19 patients (n = 35 individuals). Loess smoothing in orange. Figure 3. PMID: 32838342.

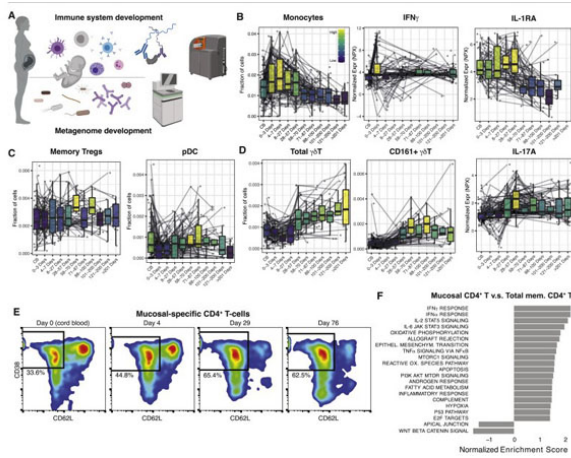

Flow Cytometry

Systems-level analysis of immune development in human newborns.

(A) Study overview.

(B) Monocyte abundance analyzed by mass cytometry and IFN γ and IL1RA measured by Olink assays in longitudinal blood samples (n = 858) from 208 individual children and binned by sampling day of life. Boxplots are colored by mean rank. CB, cord blood.

(C) Blood mass cytometry analyses of memory Tregs. pDC, plasmacytoid DC.

(D) Blood $\gamma\delta$ T cell abundance and subset of $\gamma\delta$ T cells expressing CD161 and plasma IL-17A.

(E) Representative fluorescence-activated cell sorting (FACS) plots of CD38+CD62L-CD4+ T cells sorted at postnatal days 0, 4, 29, and 76 from newborn peripheral blood mononuclear cells (PBMCs) and subjected to bulk mRNA sequencing (mRNA-seq).

(F) Gene set enrichment analysis showing the top enriched hallmark pathways in mucosal-specific versus total memory CD4+ T cells. Figure 1. PMID: 34143954.

Bottle

0.5 M EDTA pH 8.0



References

- Henrick BM et al. Bifidobacteria-mediated immune system imprinting early in life. *Cell*. (2021)
- Rodriguez L et al. Systems-Level Immunomonitoring from Acute to Recovery Phase of Severe COVID-19 *Cell Rep Med*. (2020)
- Lakshmikanth T et al. Human immune system variation during 1 year. *Cell Rep*. (2020)
- Olin A et al. Stereotypic immune system development in newborn children. *Cell*. (2018)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.